

Analysis: XHDa-re-3

Path: Home/YL Research

Based on method: Triplets Method

Created	09-Jun-2019 15:18:53	visionCATSuser
Modified	09-Jun-2019 16:58:49	visionCATSuser
Last HPTLC log	09-Jun-2019 16:58:49	Analysis modified
Explorer notes		

Track	Vial ID	Description	Volume	Position	Type
1	MeOH blank	MeOH Blank	2.0 µl	A1	Sample
2	Mixture 100	Mixture 500ng	5.0 µl	B1	Sample
3	9-THC 100	D9-THC 500ng	5.0 µl	C1	Reference
4	CBD 100	CBD 500ng	5.0 µl	D1	Reference
5	CBN 100	CBN 500ng	5.0 µl	E1	Reference
6	CBG 100	CBG 500ng	5.0 µl	F1	Reference
7	CBC 100	CBC 500ng	5.0 µl	A2	Reference
8	THCV 100	THCV 500ng	5.0 µl	B2	Reference
9	CBDV 100	CBDV 500ng	5.0 µl	C2	Reference
10	8-THC 100	D8-THC 500ng	5.0 µl	D2	Reference
11	THCA-A 100	THCA-A 500ng	5.0 µl	E2	Reference
12	CBDA 100	CBDA 500ng	5.0 µl	F2	Reference
13	CBGA 100	CBGA 500ng	5.0 µl	A3	Reference
14	Mixture 100	Mixture 500ng	5.0 µl	B1	Sample
15	MeOH blank	MeOH Blank	2.0 µl	A1	Sample

Sequence table notes

A track marked with ⚠ means: the application type is overridden in some evaluation(s).

System setup:

Software	Server User-PC, version 2.5.18072.1
ATS4	S/N:080713
Chamber	N/A
Derivatization dip	N/A
Scanner3	S/N:031025
Visualizer	S/N:230515

Chromatography

Plate layout:

Stationary phase	Merck, HPTLC plates silica gel 60 F 254
Plate format	200.0 x 100.0 mm
Application type	User
Application	Position Y: 10.0 mm, length: 8.0 mm, width: 0.0 mm
Track	First position X: 20.0 mm, distance: 11.4 mm
Solvent front position	70.0 mm
Notes	

Take image clean plate 1a - Visualizer (S/N: 230515):

XHDa-re-3

visionCATS

Quality	Enhanced
RT White	auto capture, Auto, level 85 %, Band
R 254	auto capture, Auto, level 85 %, Band
Instrument diagnostics	Valid diagnostics
Documentation step label	
Notes	

Application 1 - ATS 4 (S/N: 080713):

Spray gas	NI
Sample solvent type	Methanol
Filling speed	15 µl/s
Predosage volume	200 nl
Retraction volume	200 nl
Dosage speed	150 nl/s
Filling quality	User
Rinsing cycles / vacuum	1 / 4 s
Filling cycles / vacuum	1 / 4 s
Rinsing solvent name	Methanol
Nozzle temperature	Unheated
Rack in use	Standard
Instrument diagnostics	Valid diagnostics
Notes	

Development 1 - Chamber:

Tank	TTC 20x10
Mobile phase	
Saturation time	20 min
Use saturation pad	true
Use smartALERT	false
Volume front through	10 ml
Volume rear through	20 ml
Drying time	5 min
Drying temperature	Room temperature
Notes	

Take image developed plate 1a - Visualizer (S/N: 230515):

Quality	Enhanced
RT White	auto capture, Auto, level 85 %, Band
R 254	auto capture, Auto, level 85 %, Band
R 366	auto capture, Auto, level 85 %, Band
Instrument diagnostics	Valid diagnostics
Documentation step label	
Notes	

Scan developed plate 1b - Scanner 3 (S/N: 031025):

XHDa-re-3

visionCATS

Scanner type	Single λ
Optimization for	Resolution
Measurement mode	Absorption
Filter	n/a
Detector mode	Automatic
Scanning speed	20 mm/s
Data resolution	100 $\mu\text{m}/\text{step}$
Slit	5 x 0.2 mm, micro
Partial scan	No
Lamp	Deuterium & Tungsten
Wavelength(s)	254 nm
Instrument diagnostics	Valid diagnostics
Documentation step label	
Notes	

Derivatization 1 - dip:

Reagent name	Fast Blue B salt
Dipping speed	3
Dipping time	5 s
Reagent preparation	1g Fast Blue B salt in 200mL water
Heating	none
Notes	Air dry for 5 minutes

Take image derivatized plate 1a - Visualizer (S/N: 230515):

Quality	Enhanced
RT White	auto capture, Auto, level 85 %, Band
R 366	auto capture, Auto, level 85 %, Band
Instrument diagnostics	Valid diagnostics
Documentation step label	
Notes	

System suitability tests:

SST settings:

SST tracks	
------------	--

Data acquisition

Application 1 - ATS 4 (S/N: 080713):

Executed	09-Jun-2019 15:20:45 visionCATSuser
----------	-------------------------------------

Development 1 - Chamber:


Executed	09-Jun-2019 15:42:12 visionCATSuser
----------	-------------------------------------

Take image developed plate 1a - Visualizer (S/N: 230515):

Executed	09-Jun-2019 16:36:49 visionCATSuser
----------	-------------------------------------

XHDa-re-3
RT White


visionCATS
Developed, RemTransVis

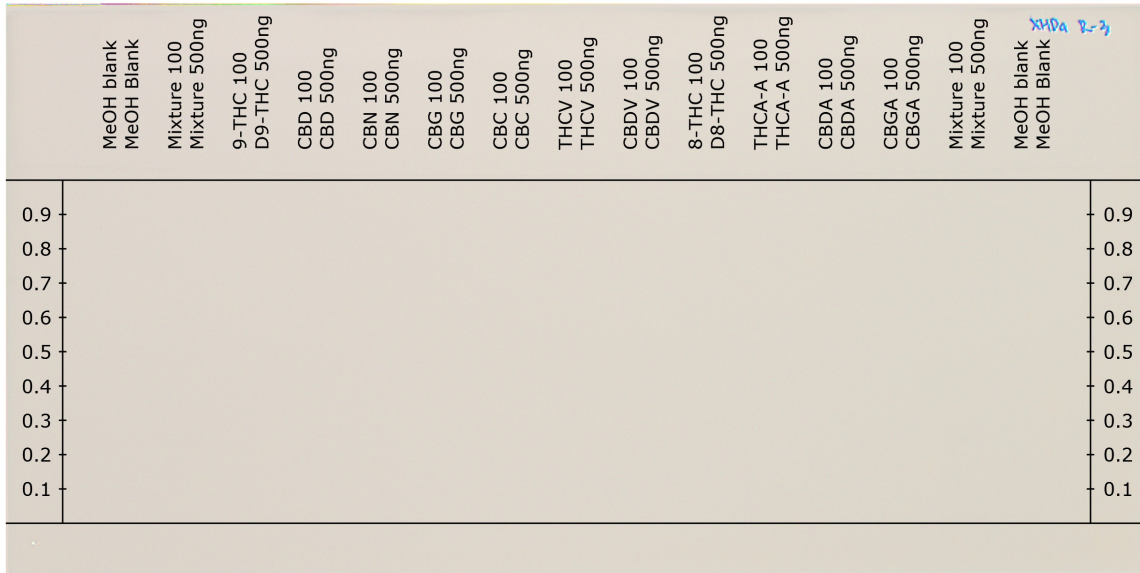

Exposure	0.085 s
Contrast	1
Normalized exposure	Disabled
Clarify	Disabled
White balance	1.00, 1.00, 1.00

R 254


Developed, Remission254

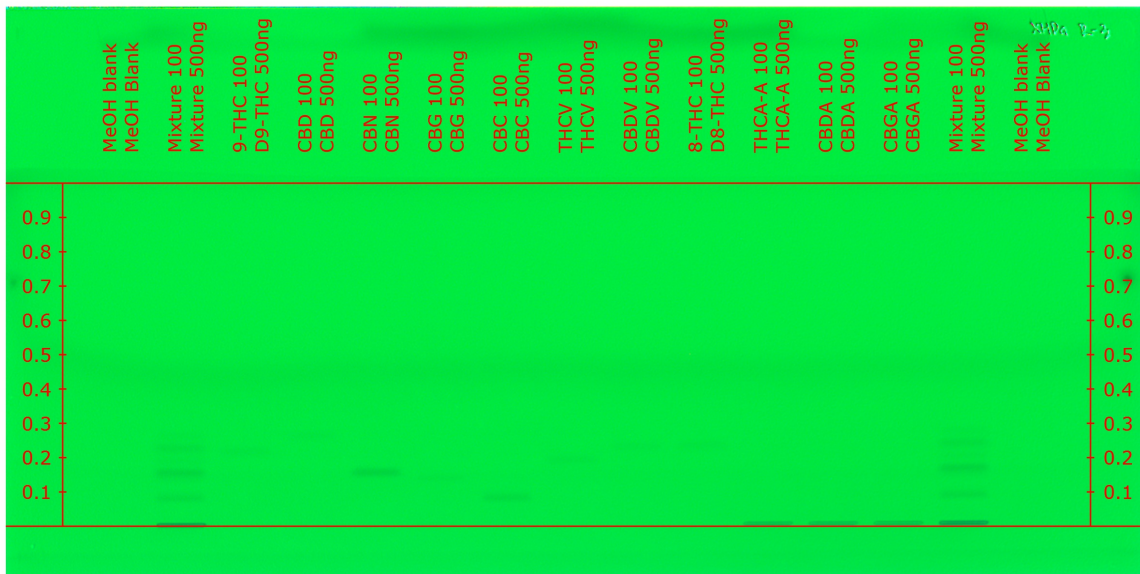

Exposure	0.264 s
Contrast	1
Normalized exposure	Disabled
Clarify	Disabled
White balance	1.00, 1.00, 1.00

XHDa-re-3
R 366

visionCATS
Developed, Remission366

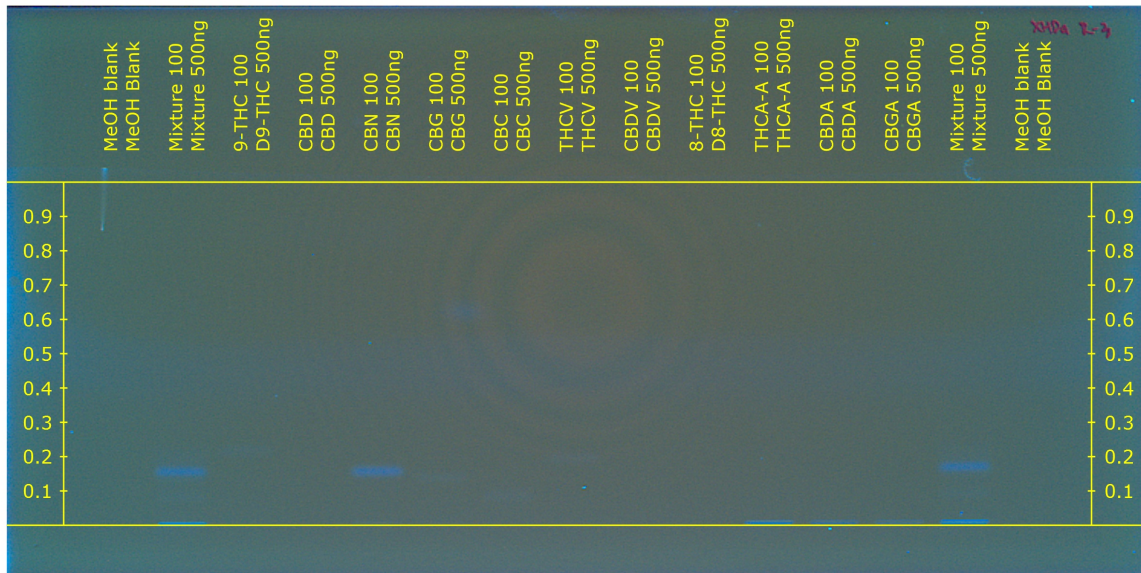

Exposure	9.999 s
Contrast	1
Normalized exposure	Disabled
Clarify	Disabled
White balance	1.00, 1.00, 1.00

Scan developed plate 1b - Scanner 3 (S/N: 031025):


Executed	09-Jun-2019 16:45:14 visionCATSuser
----------	-------------------------------------

Scan:

Wavelength	254 nm
------------	--------

Track 1:

Type	Single λ
------	------------------

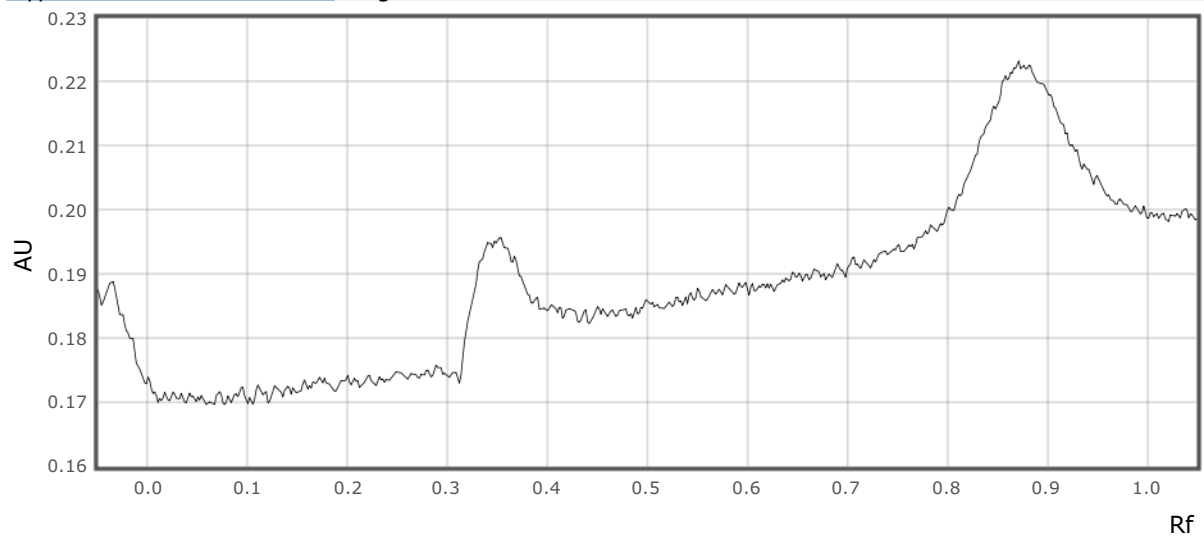


XHDa-re-3

visionCATS


Track 2:

Type Single λ

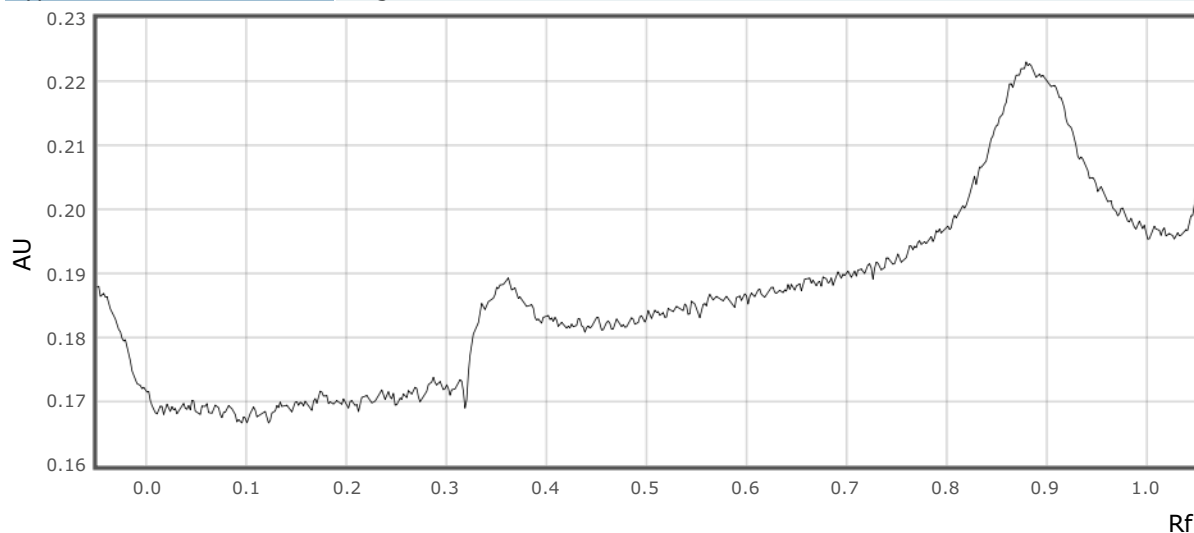


Track 3:

Type Single λ

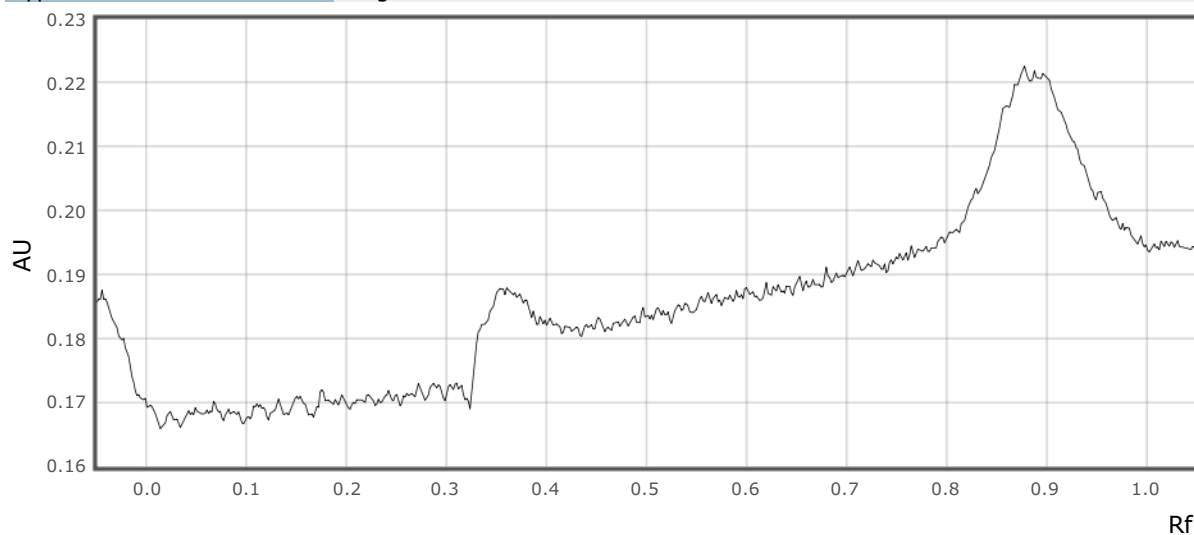


Track 4:

Type Single λ

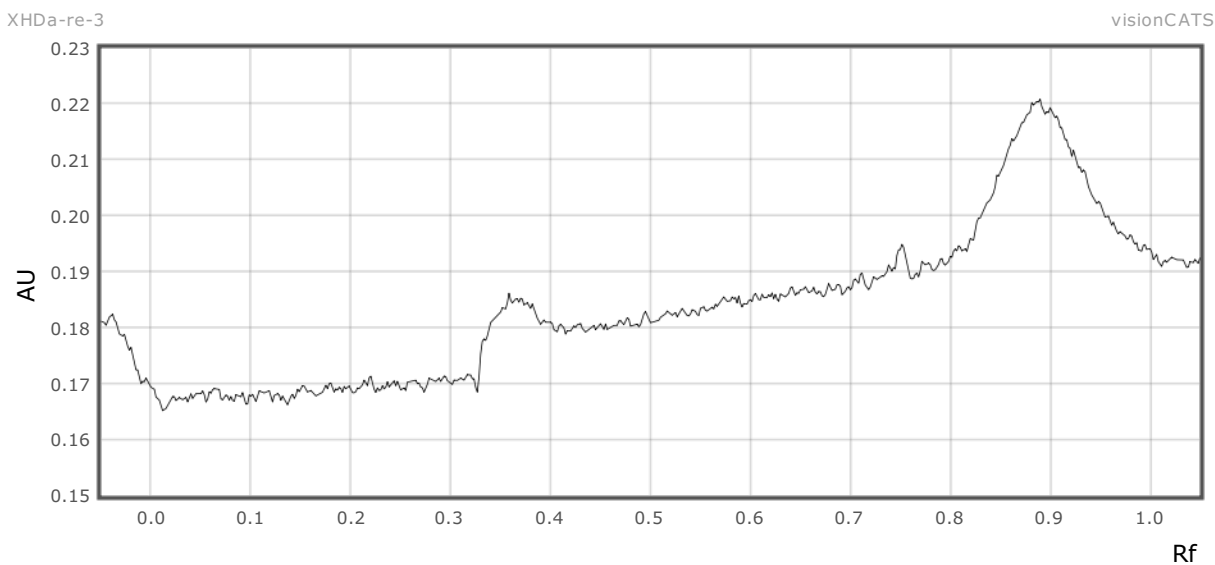


Track 5:

Type Single λ

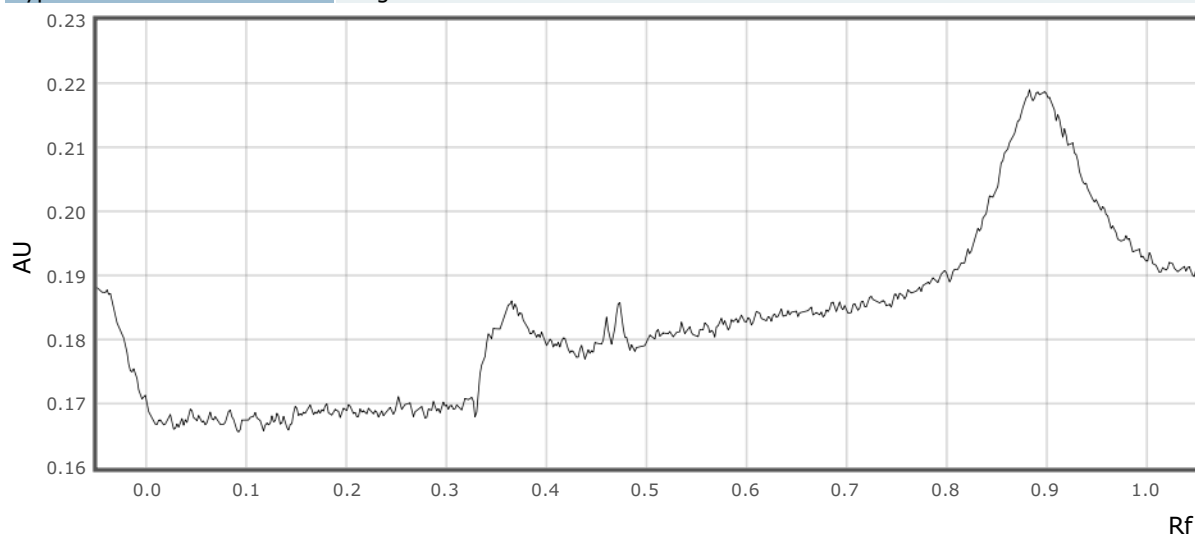


Track 6:

Type Single λ

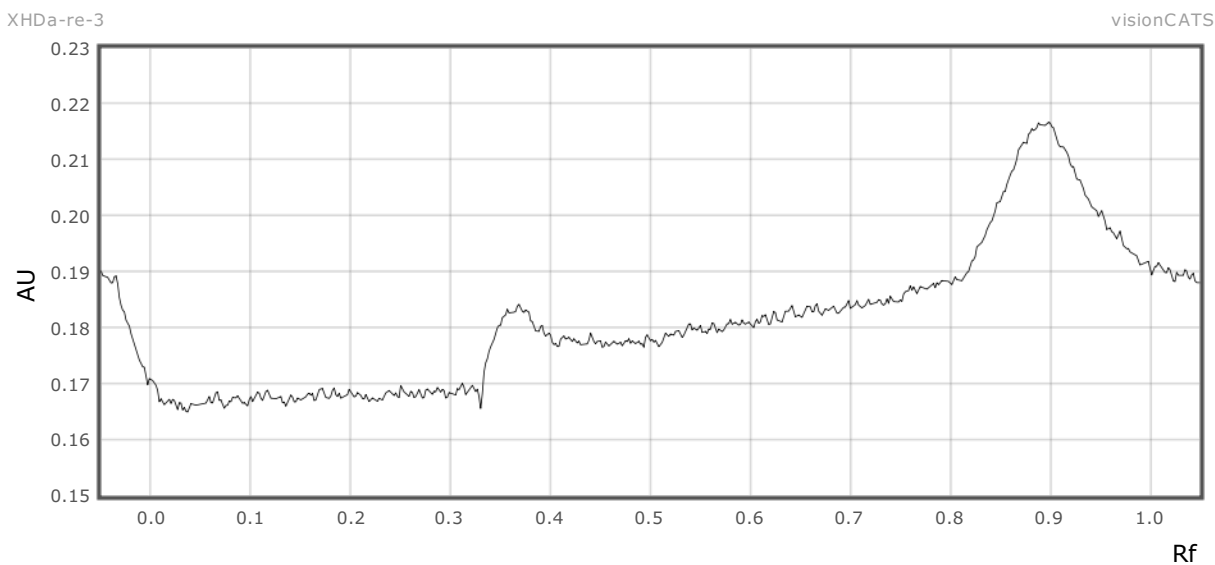


Track 7:

Type Single λ

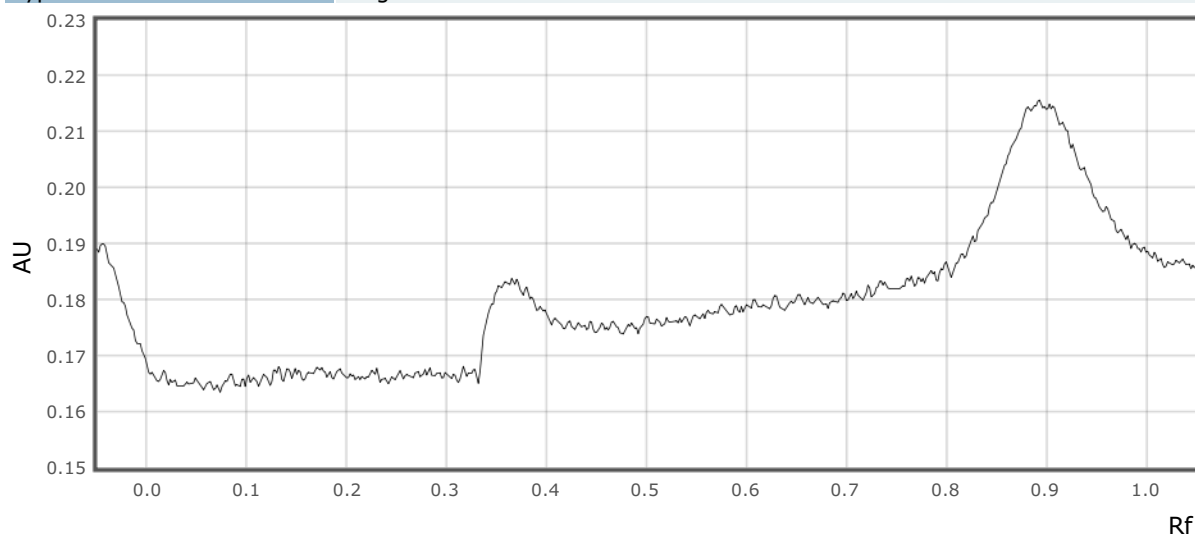


Track 8:

Type Single λ

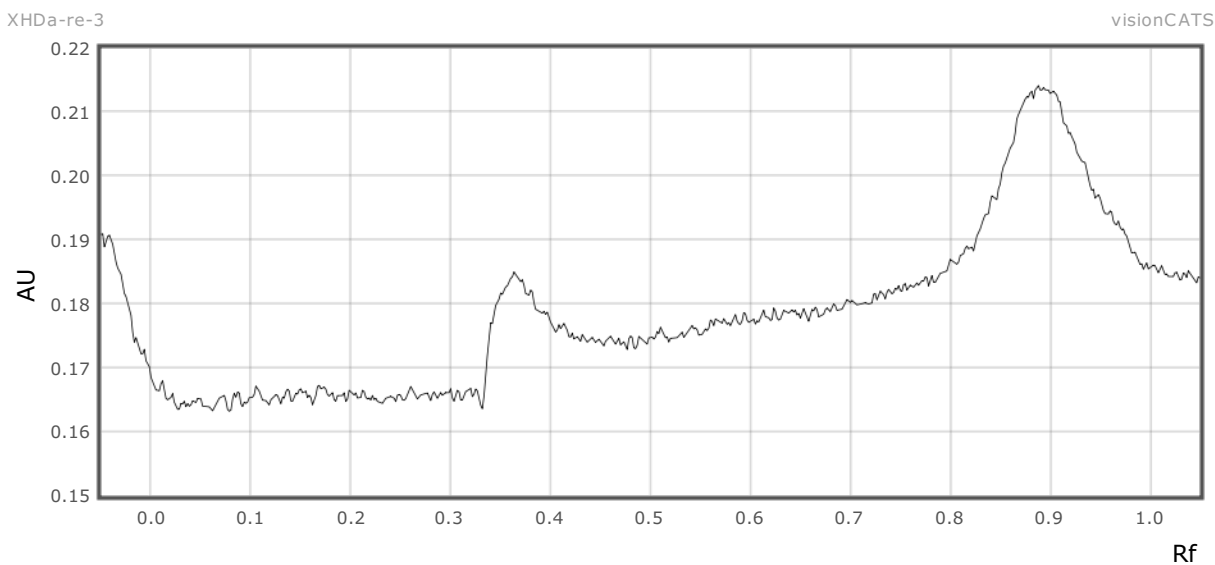


Track 9:

Type Single λ

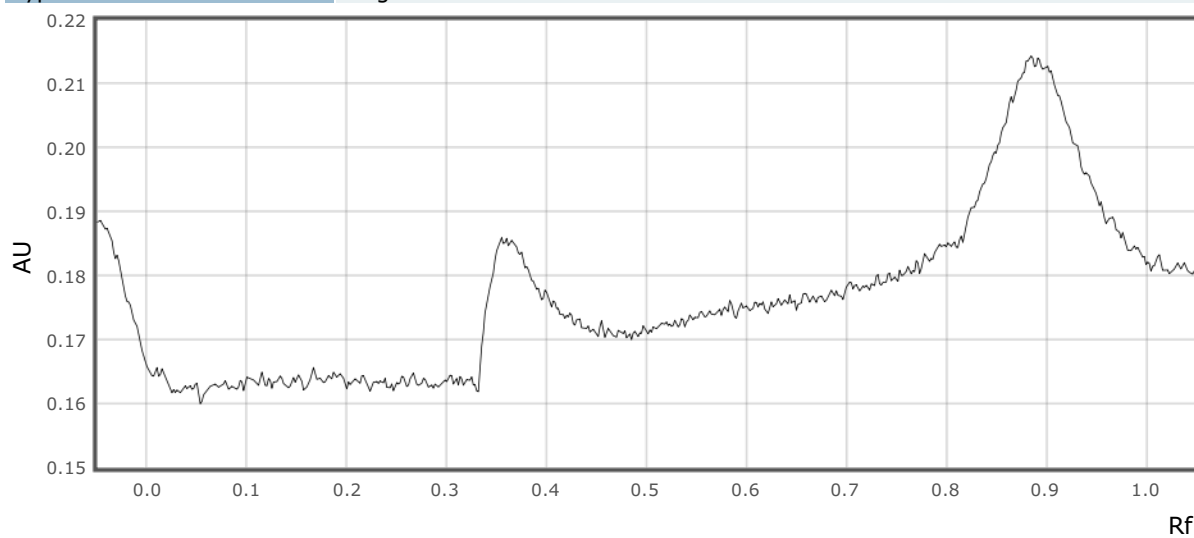


Track 10:

Type Single λ

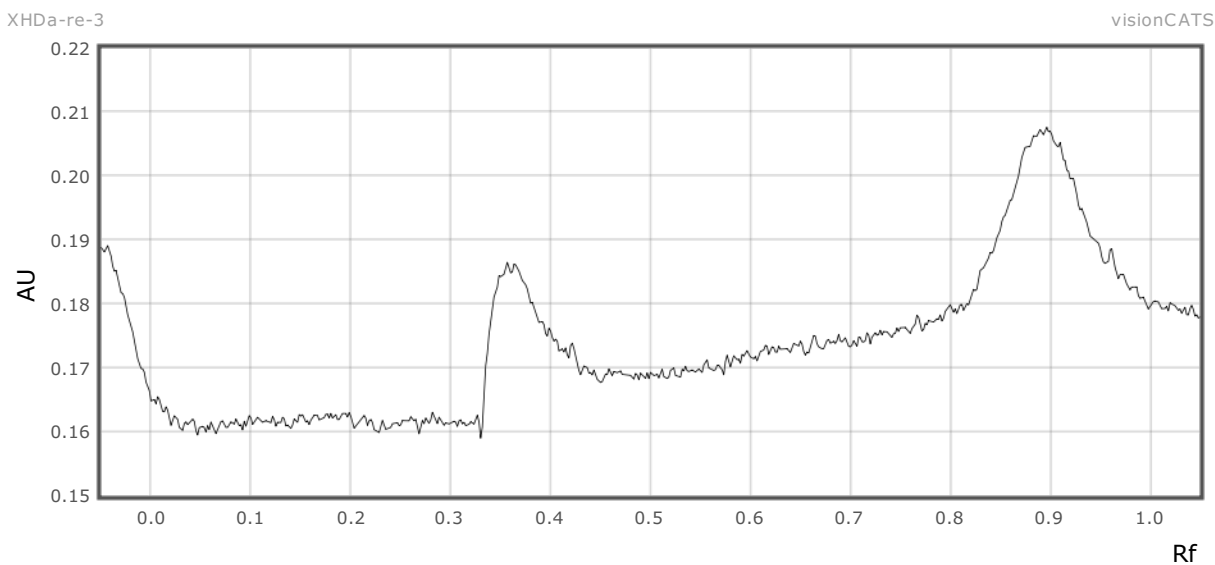

Track 11:

Type Single λ

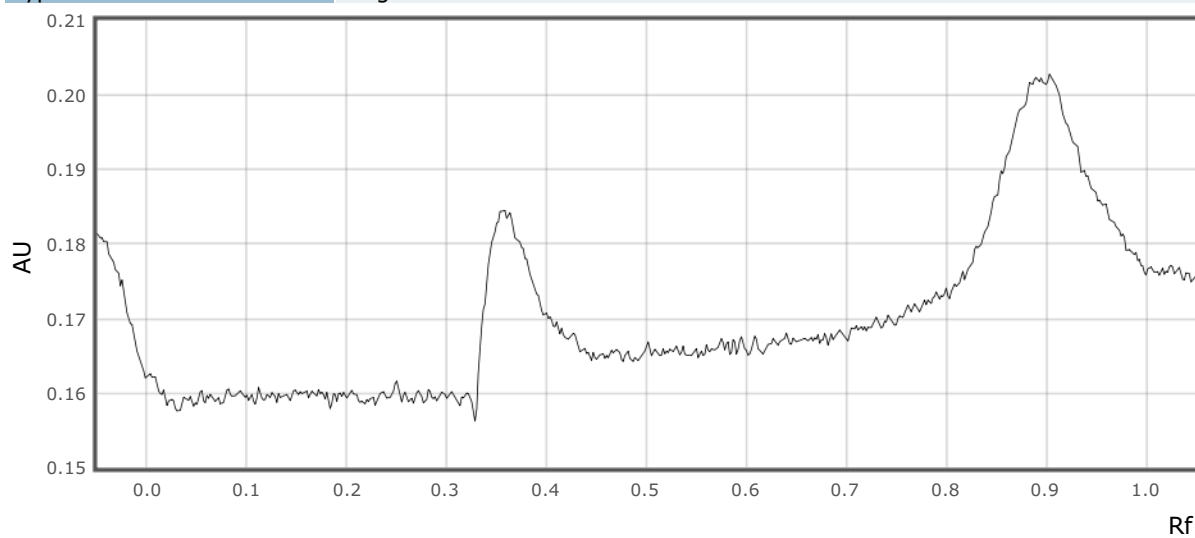


Track 12:

Type Single λ


XHDa-re-3

visionCATS

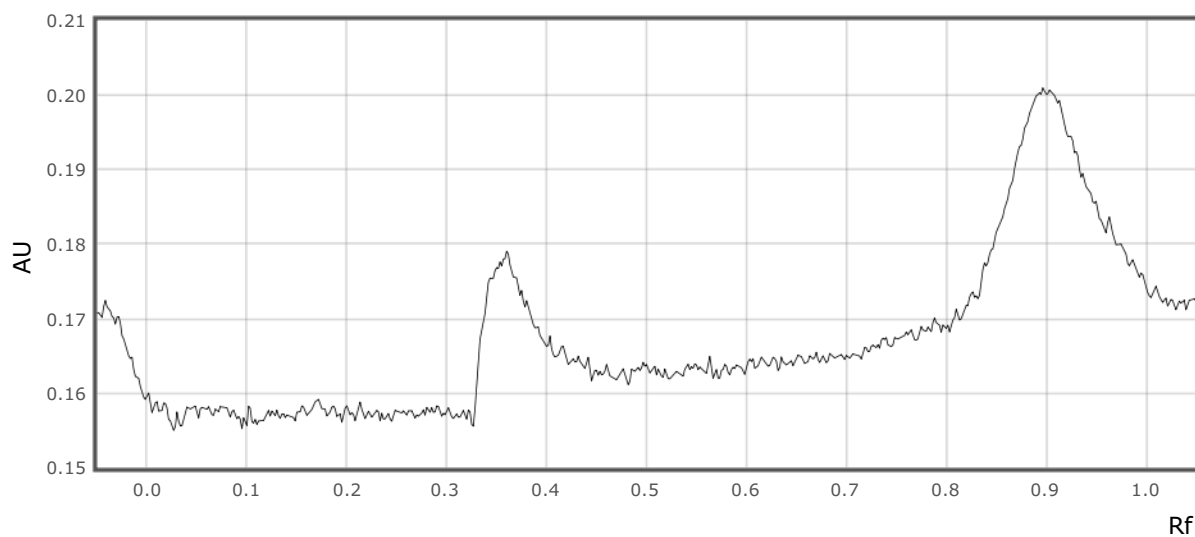


Track 13:

Type Single λ

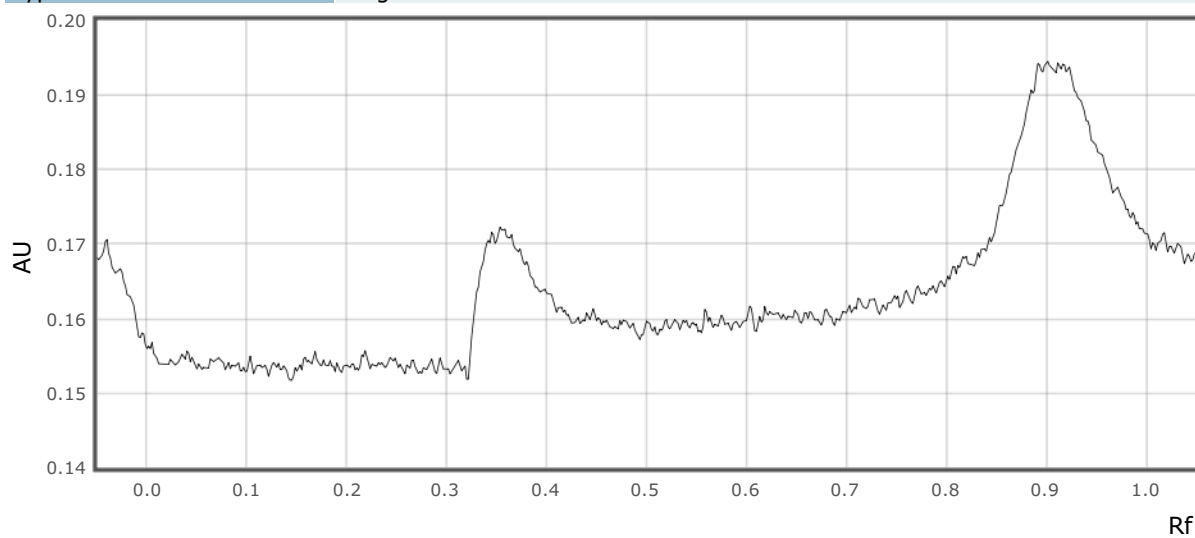


Track 14:

Type Single λ

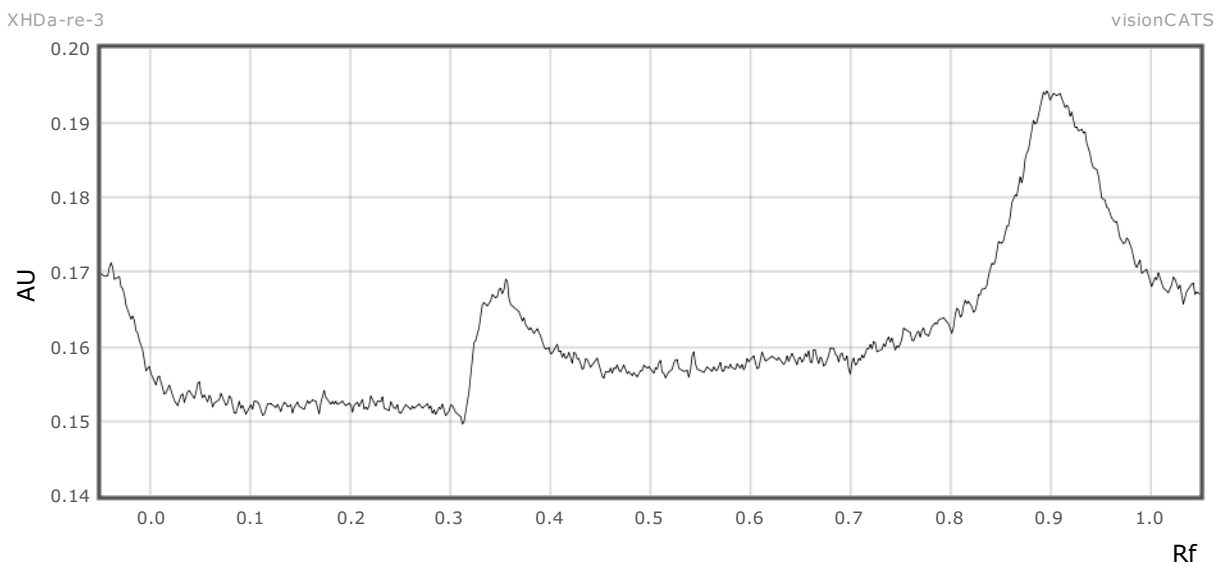

Track 15:

Type Single λ

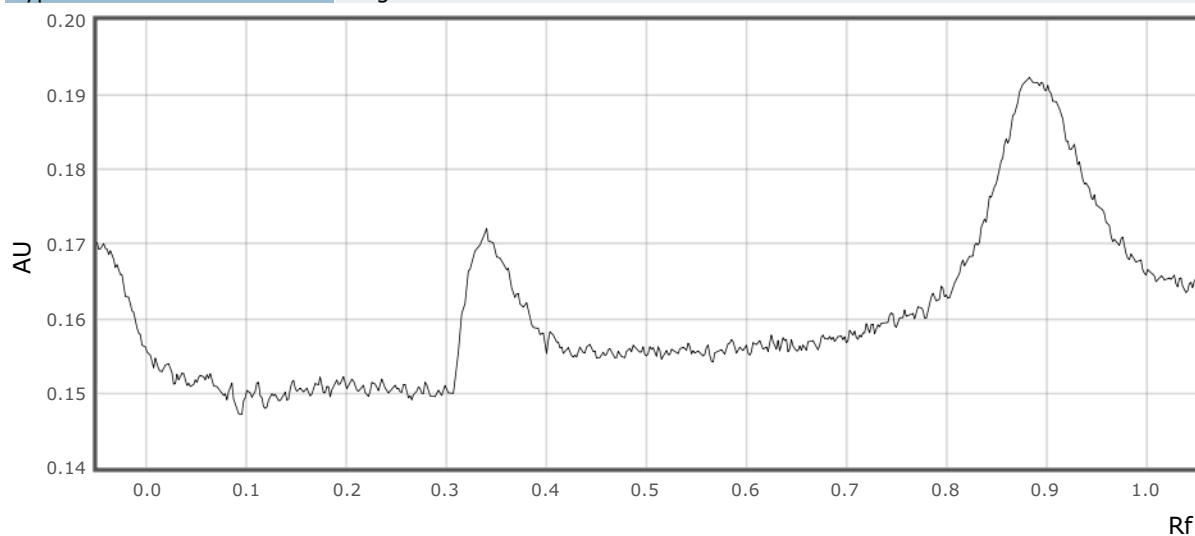

Derivatization 1 - dip:


Executed 09-Jun-2019 16:52:59 visionCATSuser

Take image derivatized plate 1a - Visualizer (S/N: 230515):

Executed 09-Jun-2019 16:53:02 visionCATSuser

XHDa-re-3
RT White


visionCATS
Derivatized, RemTransVis

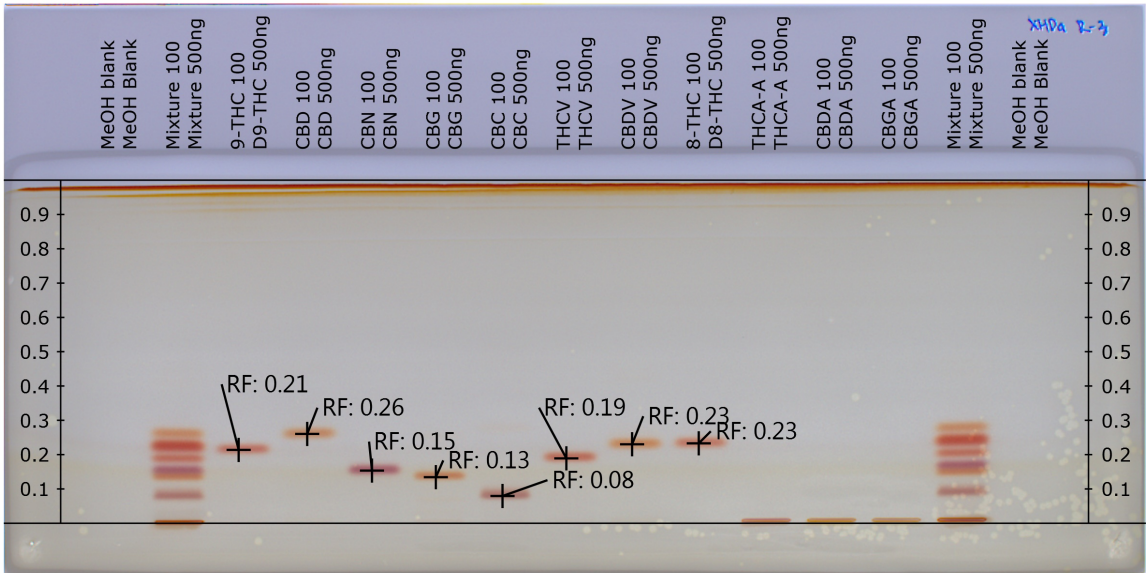

Exposure	0.046 s
Contrast	1
Normalized exposure	Disabled
Clarify	Disabled
White balance	1.20, 1.12, 0.79

R 366

Derivatized, Remission366

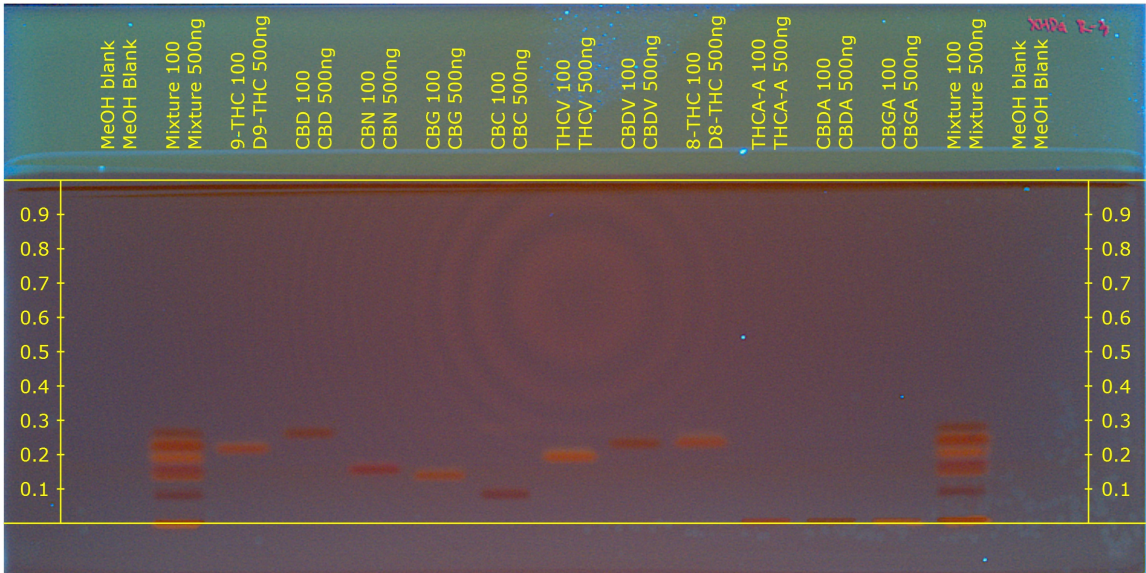

Exposure	9.999 s
Contrast	1
Normalized exposure	Disabled
Clarify	Disabled
White balance	1.00, 1.00, 1.00

Evaluation 1 :

XHDa-re-3

visionCATS

Validated	false
Step	Take image derivatized plate 1a
Concentration unit type	Mass / volume
Notes	

Definition:

References:

9-THC 100

Substance Name	Concentration	Purity
9-THC	100.000 µg/ml	100.00 %

CBD 100

Substance Name	Concentration	Purity
CBD	100.000 µg/ml	100.00 %

CBN 100

Substance Name	Concentration	Purity
CBN	100.000 µg/ml	100.00 %

CBG 100

Substance Name	Concentration	Purity
CBG	100.000 µg/ml	100.00 %

CBC 100

Substance Name	Concentration	Purity
CBC	100.000 µg/ml	100.00 %

THCV 100

Substance Name	Concentration	Purity
THCV	100.000 µg/ml	100.00 %

CBDV 100

Substance Name	Concentration	Purity
CBDV	100.000 µg/ml	100.00 %

8-THC 100

Substance Name	Concentration	Purity
8-THC	100.000 µg/ml	100.00 %

THCA-A 100

Substance Name	Concentration	Purity
THCA-A	100.000 µg/ml	100.00 %

CBDA 100

Substance Name	Concentration	Purity
CBDA	100.000 µg/ml	100.00 %

CBGA 100

Substance Name	Concentration	Purity
CBGA	100.000 µg/ml	100.00 %

XHDa-re-3


visionCATS

Samples:				
Vial ID	Amount	Volume solution	Reference amount	Related to
MeOH blank		0.00 ml		
Mixture 100		0.00 ml		

Integration parameters:	
Bounds	[0.000,1.000]
Smoothing	Savitzky-Golay of order 3 and window 7
Baseline correction	Lowest slope with noise 0.05
Profile subtraction	Profile subtraction from track 1
Peaks detection	Gauss (legacy) with sensitivity 0.1, separation 1 and threshold 0.1

Scan:	
Wavelength	RT White

Track 1:	
Type	Sample
Vial ID	MeOH blank
Description	MeOH Blank
Volume	2.0 µl

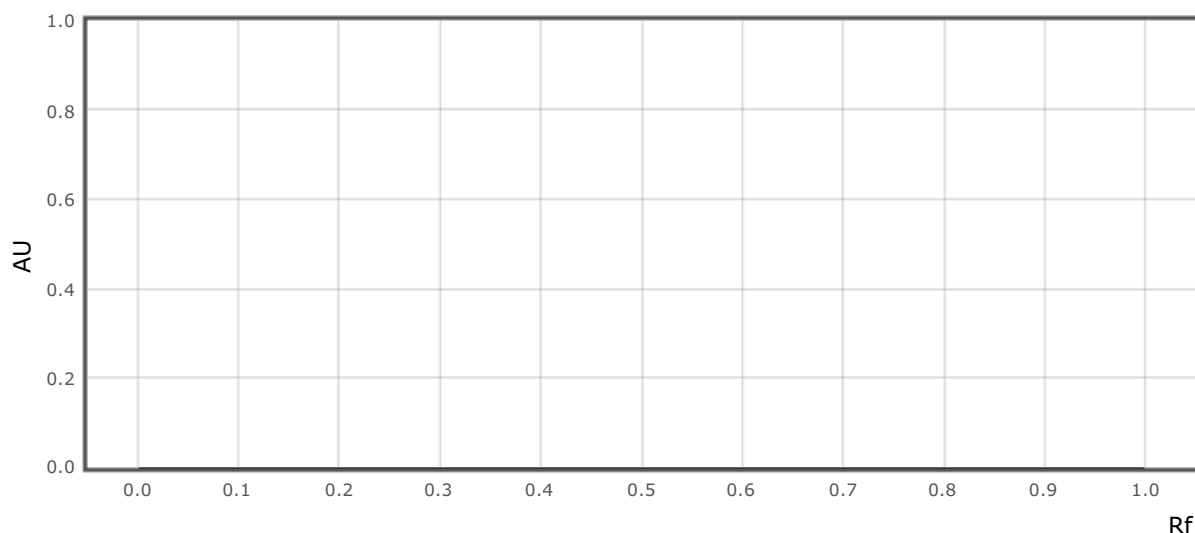


Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		

Track 2:	
Type	Sample
Vial ID	Mixture 100
Description	Mixture 500ng
Volume	5.0 µl

XHDa-re-3


visionCATS

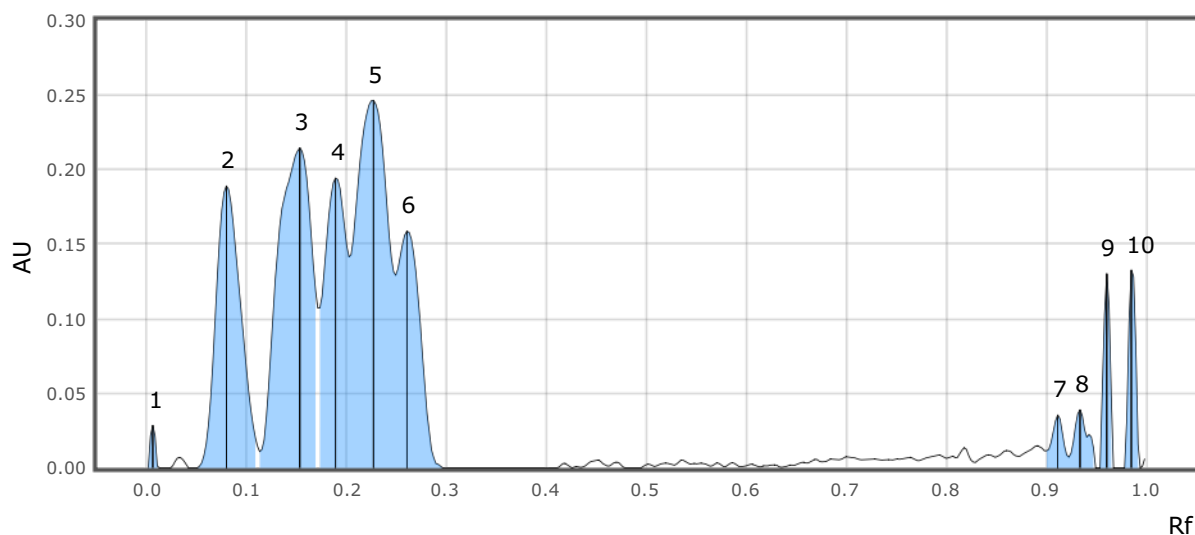

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.001	0.0000	0.006	0.0284	2.08	0.012	0.0000	0.00016	0.46	No	
2	0.050	0.0000	0.079	0.1888	13.81	0.111	0.0144	0.00539	15.36	No	
3	0.113	0.0108	0.153	0.2143	15.67	0.171	0.1072	0.00840	23.93	No	
4	0.173	0.1071	0.189	0.1941	14.20	0.202	0.1413	0.00473	13.48	No	
5	0.202	0.1413	0.227	0.2462	18.01	0.249	0.1290	0.00924	26.31	No	
6	0.249	0.1290	0.260	0.1585	11.59	0.298	0.0000	0.00405	11.53	No	
7	0.901	0.0114	0.912	0.0354	2.59	0.923	0.0071	0.00047	1.33	No	
8	0.923	0.0071	0.934	0.0389	2.84	0.950	0.0000	0.00060	1.72	No	
9	0.954	0.0000	0.961	0.1300	9.51	0.967	0.0000	0.00098	2.79	No	
10	0.979	0.0000	0.985	0.1325	9.69	0.994	0.0000	0.00109	3.10	No	

Track 3:

Type	Reference
Vial ID	9-THC 100
Description	D9-THC 500ng
Volume	5.0 µl

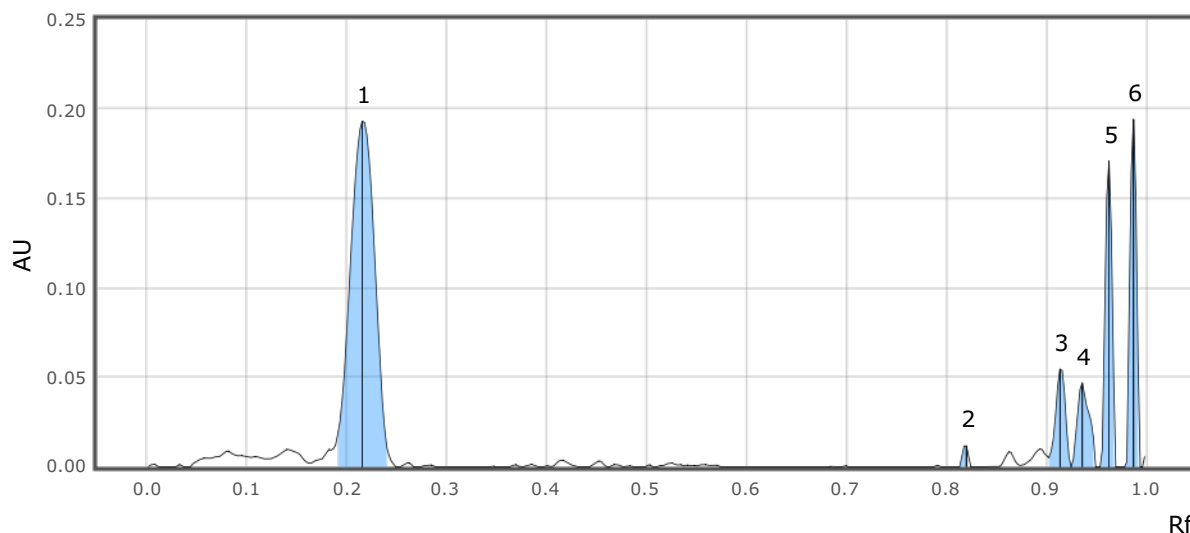


XHDa-re-3

visionCATS

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.191	0.0182	0.216	0.1929	28.76	0.240	0.0107	0.00541	56.77	Yes	9-THC
2	0.813	0.0000	0.820	0.0118	1.75	0.825	0.0000	0.00008	0.80	No	
3	0.903	0.0051	0.914	0.0546	8.14	0.925	0.0000	0.00061	6.44	No	
4	0.925	0.0000	0.936	0.0469	6.99	0.950	0.0000	0.00064	6.71	No	
5	0.954	0.0000	0.963	0.1706	25.43	0.970	0.0000	0.00132	13.81	No	
6	0.979	0.0000	0.988	0.1940	28.92	0.994	0.0000	0.00147	15.46	No	

Track 4:


Type	Reference
Vial ID	CBD 100
Description	CBD 500ng
Volume	5.0 µl

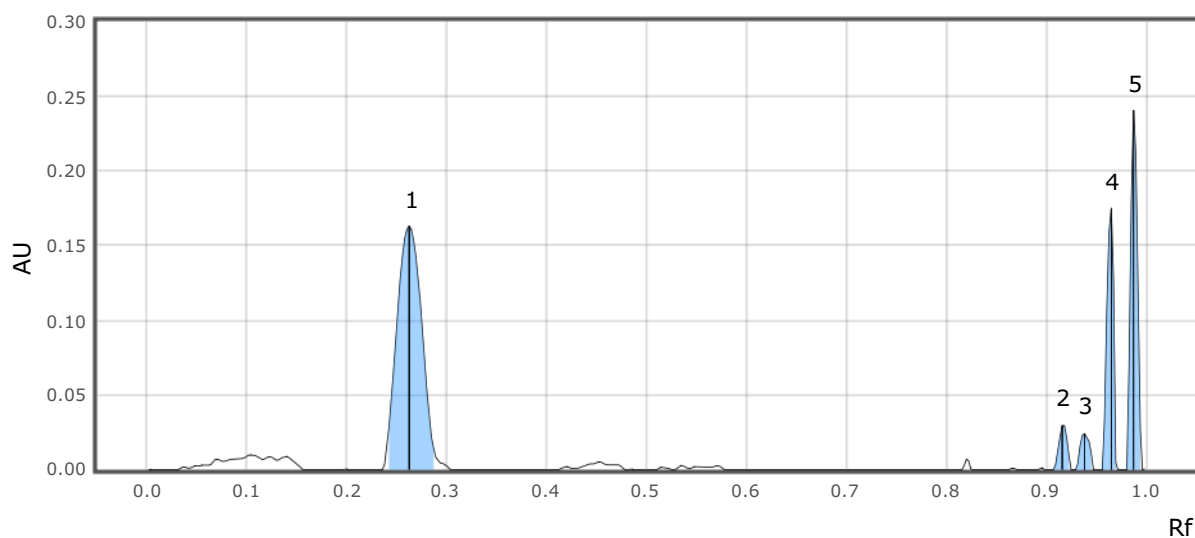

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.239	0.0156	0.262	0.1629	25.79	0.287	0.0141	0.00462	55.04	Yes	CBD
2	0.907	0.0000	0.916	0.0295	4.68	0.925	0.0000	0.00028	3.37	No	
3	0.930	0.0000	0.938	0.0242	3.83	0.947	0.0000	0.00026	3.14	No	
4	0.956	0.0000	0.965	0.1746	27.65	0.972	0.0000	0.00135	16.07	No	
5	0.981	0.0000	0.988	0.2404	38.06	0.996	0.0000	0.00188	22.39	No	

Track 5:


Type	Reference
Vial ID	CBN 100
Description	CBN 500ng
Volume	5.0 µl

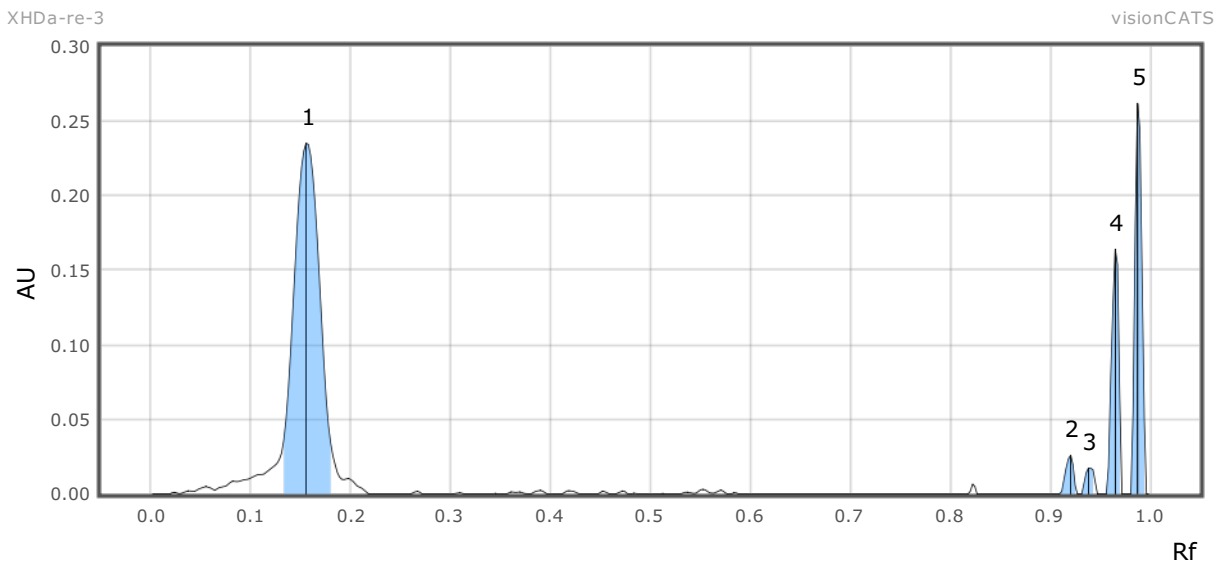

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.133	0.0361	0.155	0.2349	33.41	0.182	0.0255	0.00692	64.23	Yes	CBN
2	0.909	0.0000	0.921	0.0258	3.67	0.927	0.0000	0.00023	2.13	No	
3	0.932	0.0000	0.938	0.0172	2.45	0.947	0.0000	0.00018	1.65	No	
4	0.956	0.0000	0.965	0.1637	23.29	0.972	0.0000	0.00135	12.53	No	
5	0.981	0.0000	0.988	0.2615	37.18	0.996	0.0000	0.00210	19.47	No	

Track 6:

Type	Reference
Vial ID	CBG 100
Description	CBG 500ng
Volume	5.0 µl

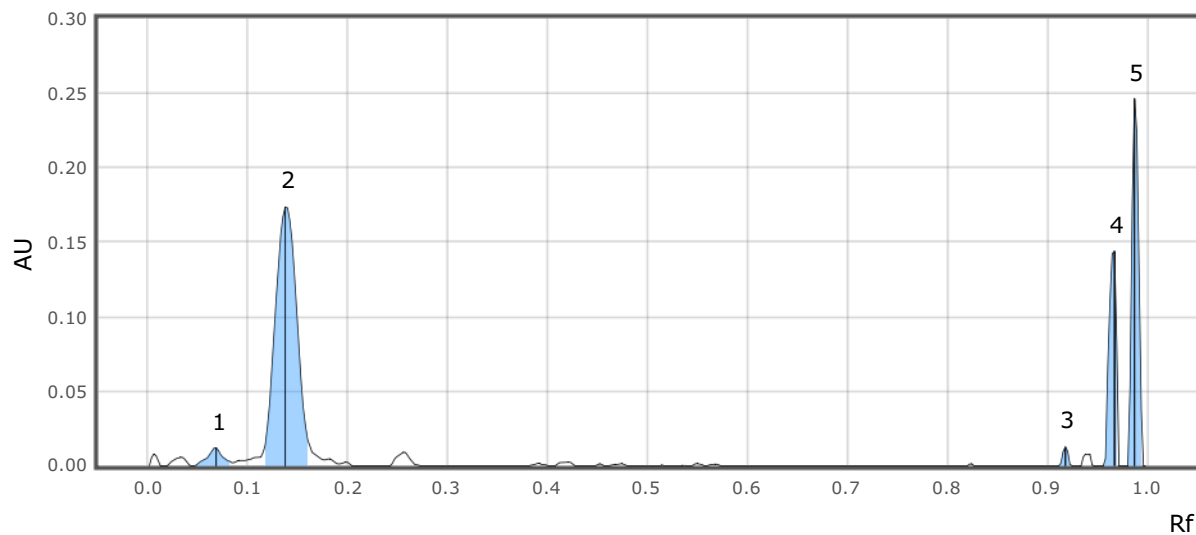


XHDa-re-3

visionCATS

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.046	0.0000	0.068	0.0122	2.07	0.084	0.0021	0.00021	2.73	No	
2	0.116	0.0089	0.137	0.1734	29.49	0.161	0.0131	0.00432	55.65	Yes	CBG
3	0.912	0.0000	0.918	0.0126	2.15	0.925	0.0000	0.00008	1.00	No	
4	0.956	0.0000	0.967	0.1438	24.46	0.972	0.0000	0.00123	15.80	No	
5	0.981	0.0000	0.988	0.2459	41.83	0.996	0.0000	0.00192	24.80	No	

Track 7:


Type	Reference
Vial ID	CBC 100
Description	CBC 500ng
Volume	5.0 µl

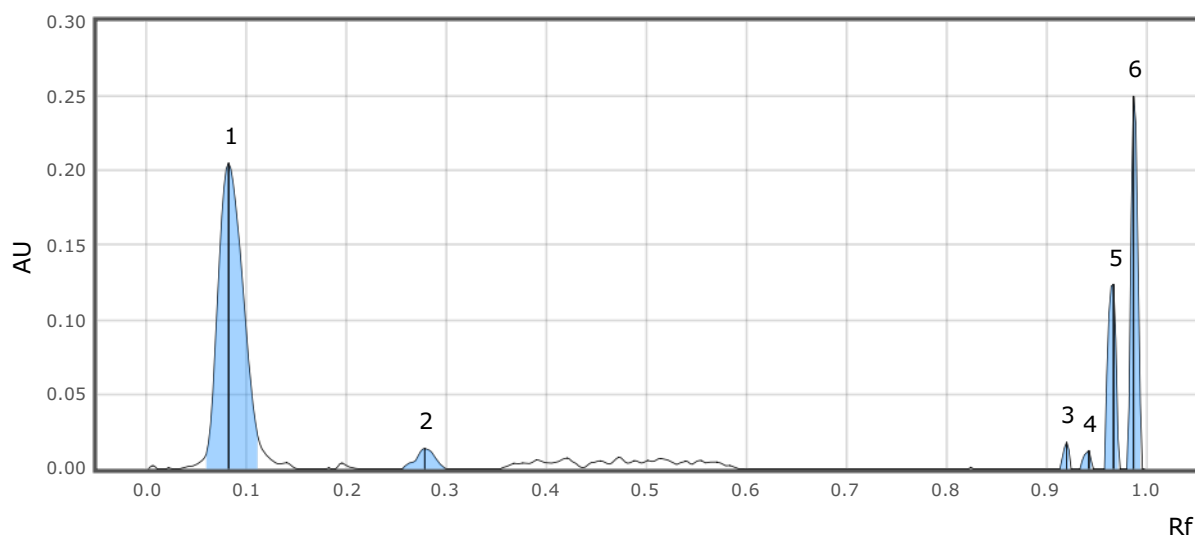

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.059	0.0096	0.082	0.2048	32.94	0.112	0.0170	0.00592	62.05	Yes	CBC
2	0.256	0.0000	0.278	0.0137	2.20	0.300	0.0000	0.00030	3.10	No	
3	0.914	0.0000	0.921	0.0178	2.86	0.925	0.0000	0.00011	1.20	No	
4	0.934	0.0000	0.943	0.0123	1.98	0.947	0.0000	0.00010	1.02	No	
5	0.959	0.0000	0.967	0.1237	19.88	0.974	0.0000	0.00116	12.17	No	
6	0.981	0.0000	0.988	0.2496	40.14	0.996	0.0000	0.00195	20.46	No	

Track 8:


Type	Reference
Vial ID	THCV 100
Description	THCV 500ng
Volume	5.0 µl

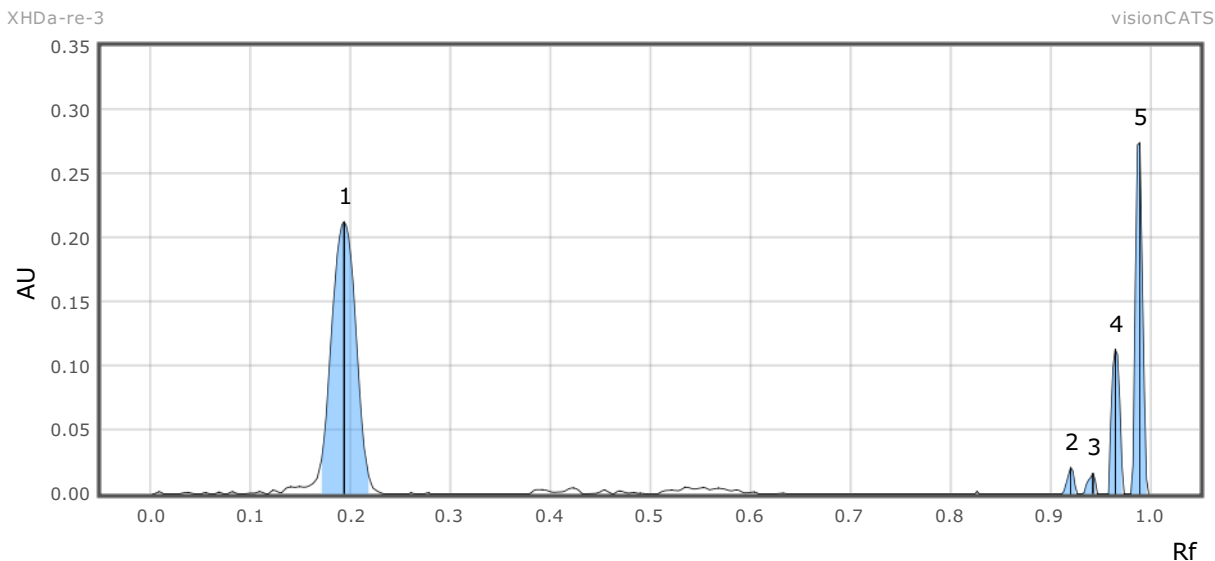

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.168	0.0192	0.193	0.2126	33.40	0.218	0.0144	0.00596	62.36	Yes	THCV
2	0.912	0.0000	0.921	0.0206	3.24	0.927	0.0000	0.00016	1.63	No	
3	0.932	0.0000	0.943	0.0160	2.51	0.947	0.0000	0.00013	1.35	No	
4	0.959	0.0000	0.965	0.1129	17.73	0.974	0.0000	0.00105	11.03	No	
5	0.981	0.0000	0.990	0.2745	43.12	0.999	0.0000	0.00226	23.63	No	

Track 9:

Type	Reference
Vial ID	CBDV 100
Description	CBDV 500ng
Volume	5.0 µl

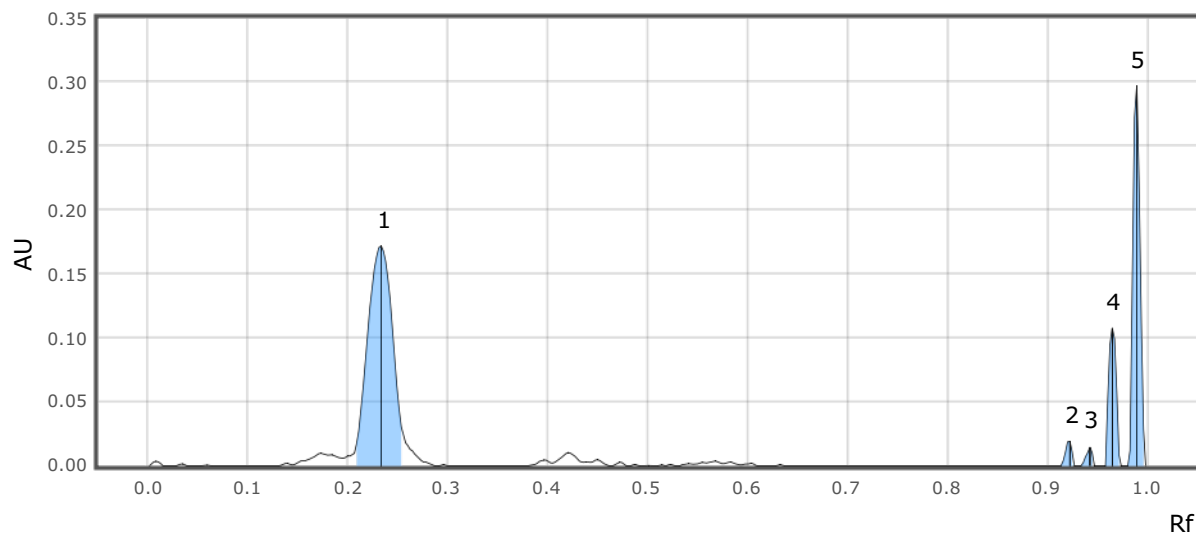


XHDa-re-3

visionCATS

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.208	0.0172	0.233	0.1719	28.18	0.258	0.0179	0.00496	57.88	Yes	CBDV
2	0.914	0.0000	0.923	0.0191	3.13	0.927	0.0000	0.00014	1.63	No	
3	0.934	0.0000	0.943	0.0145	2.37	0.947	0.0000	0.00010	1.15	No	
4	0.959	0.0000	0.965	0.1076	17.64	0.974	0.0000	0.00094	10.94	No	
5	0.981	0.0000	0.990	0.2970	48.69	0.999	0.0000	0.00243	28.40	No	

Track 10:


Type	Reference
Vial ID	8-THC 100
Description	D8-THC 500ng
Volume	5.0 µl

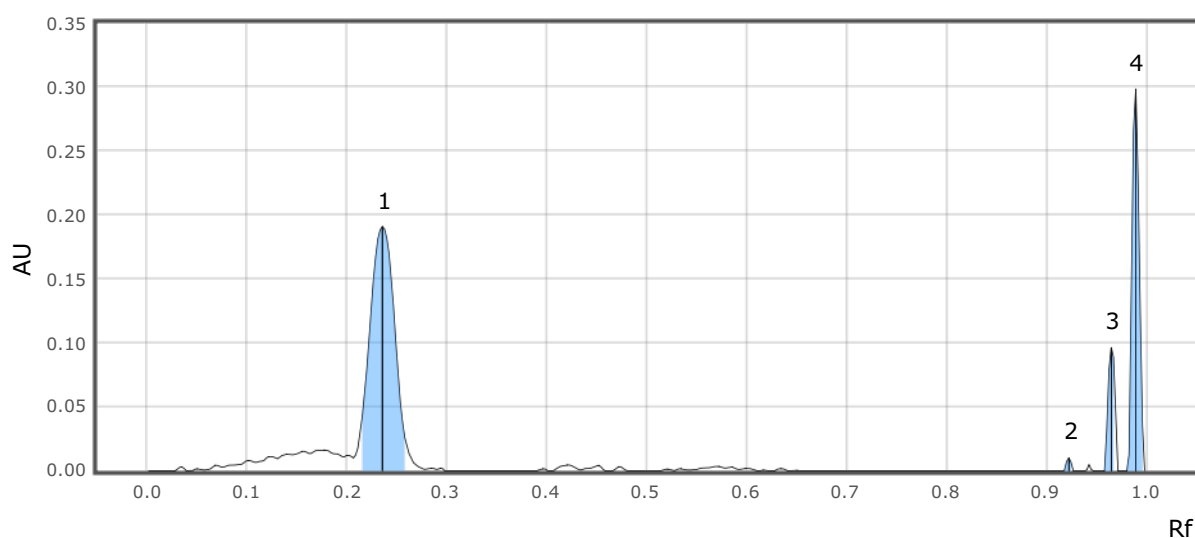

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.213	0.0296	0.236	0.1908	32.04	0.261	0.0198	0.00541	62.31	Yes	8-THC
2	0.916	0.0000	0.923	0.0101	1.69	0.927	0.0000	0.00005	0.62	No	
3	0.959	0.0000	0.965	0.0963	16.17	0.972	0.0000	0.00076	8.70	No	
4	0.981	0.0000	0.990	0.2985	50.10	0.999	0.0000	0.00246	28.37	No	

Track 11:


Type	Reference
Vial ID	THCA-A 100
Description	THCA-A 500ng
Volume	5.0 µl

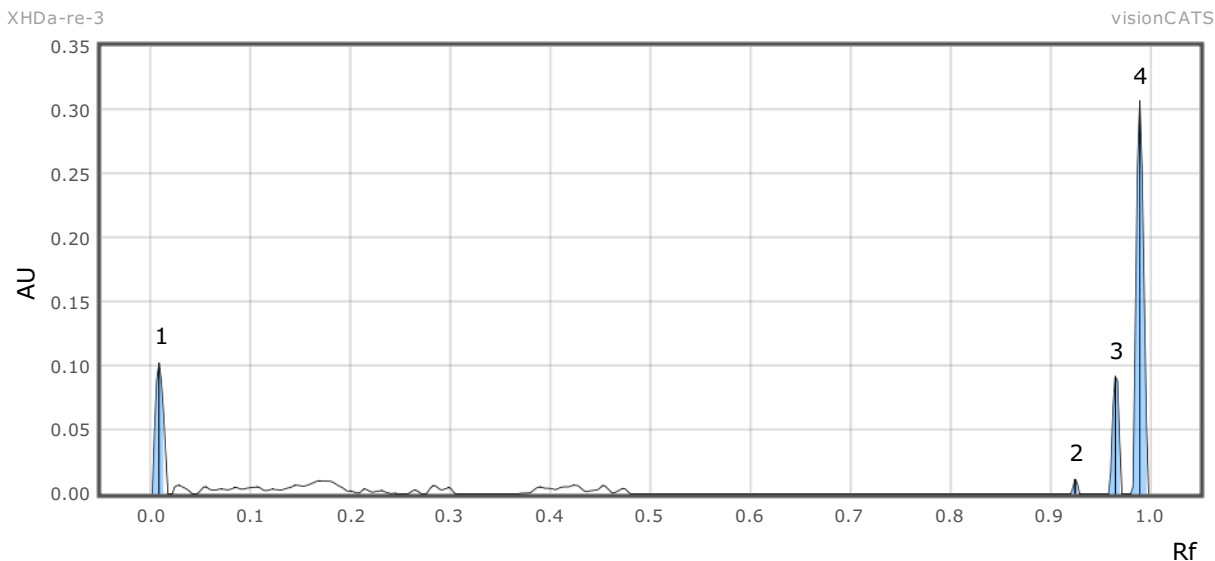

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.001	0.0000	0.008	0.1024	19.96	0.017	0.0000	0.00095	22.18	No	THCA-A
2	0.921	0.0000	0.925	0.0115	2.25	0.930	0.0000	0.00005	1.26	No	
3	0.959	0.0000	0.965	0.0917	17.88	0.972	0.0000	0.00068	15.89	No	
4	0.981	0.0000	0.990	0.3073	59.91	0.999	0.0000	0.00259	60.67	No	

Track 12:

Type	Reference
Vial ID	CBDA 100
Description	CBDA 500ng
Volume	5.0 µl

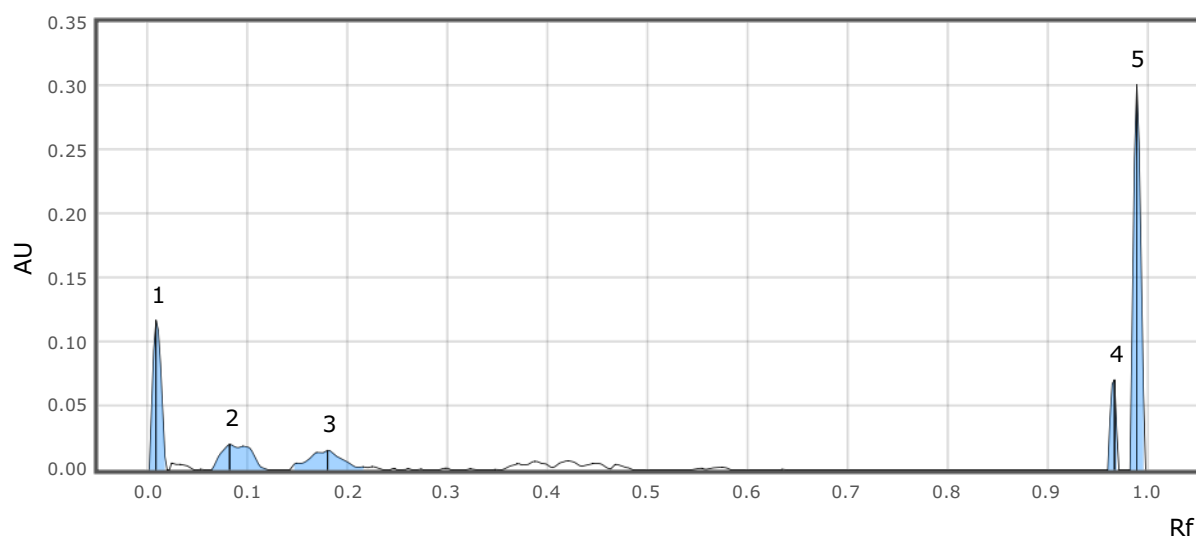


XHDa-re-3

visionCATS

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.001	0.0000	0.008	0.1169	22.33	0.019	0.0000	0.00114	21.09	No	CBDA
2	0.064	0.0000	0.082	0.0200	3.81	0.122	0.0000	0.00069	12.69	No	
3	0.146	0.0044	0.180	0.0152	2.91	0.220	0.0017	0.00059	10.92	No	
4	0.961	0.0000	0.967	0.0704	13.44	0.972	0.0000	0.00046	8.59	No	
5	0.983	0.0000	0.990	0.3010	57.51	0.999	0.0000	0.00252	46.72	No	

Track 13:

Type	Reference
Vial ID	CBGA 100
Description	CBGA 500ng
Volume	5.0 µl

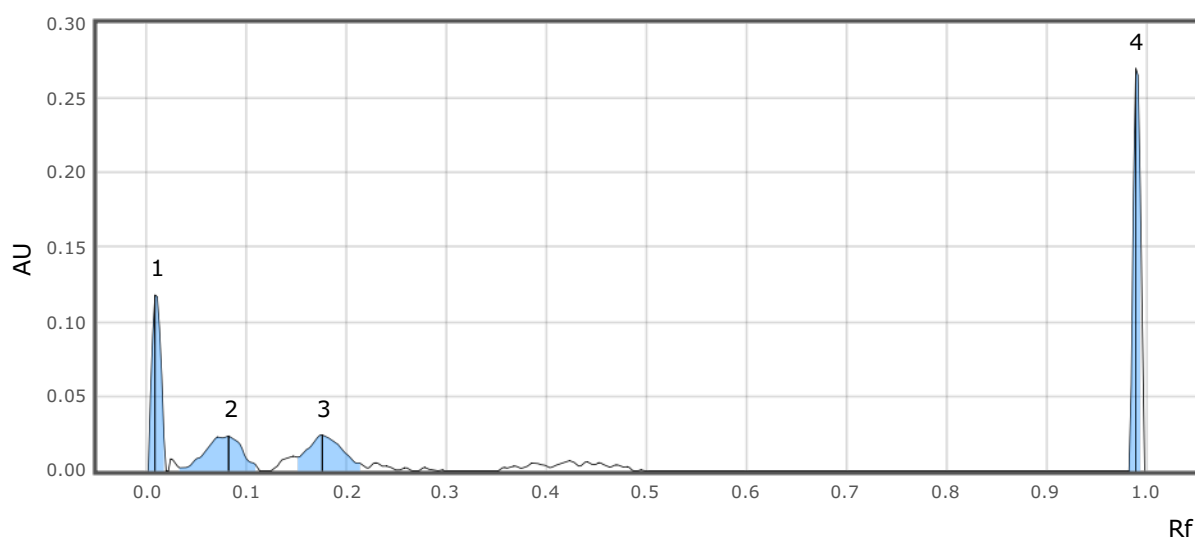


Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.001	0.0000	0.008	0.1178	27.09	0.019	0.0000	0.00124	22.29	No	CBGA
2	0.033	0.0020	0.082	0.0234	5.37	0.113	0.0000	0.00100	17.96	No	
3	0.151	0.0092	0.175	0.0240	5.53	0.220	0.0019	0.00098	17.56	No	
4	0.983	0.0000	0.990	0.2697	62.01	0.999	0.0000	0.00235	42.18	No	

Track 14:

Type	Sample
Vial ID	Mixture 100
Description	Mixture 500ng
Volume	5.0 µl

XHDa-re-3


visionCATS

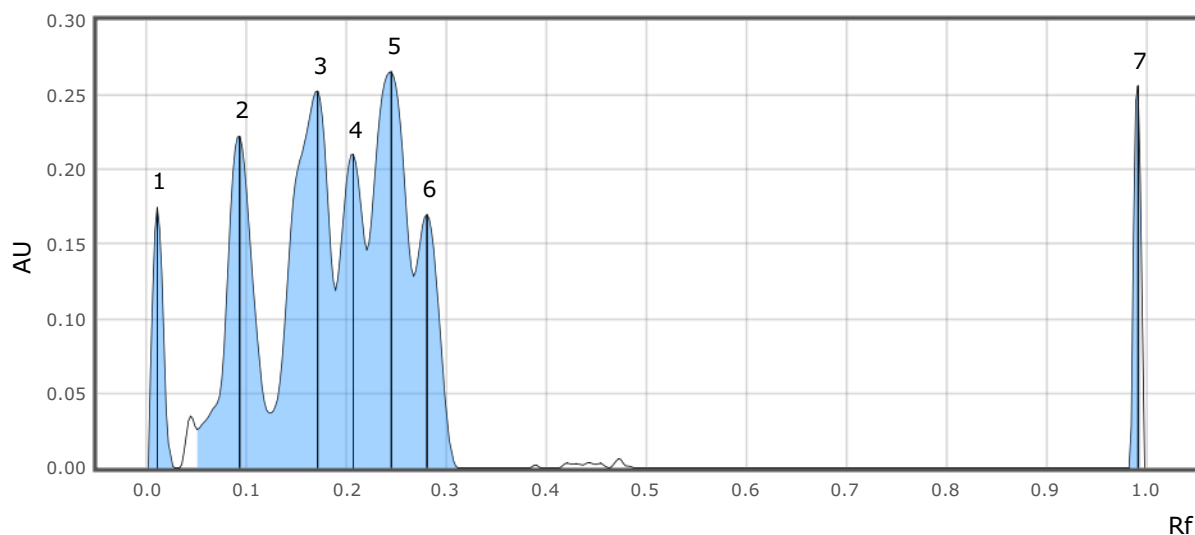

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.001	0.0000	0.010	0.1745	11.26	0.028	0.0000	0.00214	5.07	No	
2	0.050	0.0256	0.093	0.2220	14.32	0.124	0.0367	0.00730	17.34	No	
3	0.124	0.0367	0.171	0.2524	16.29	0.189	0.1187	0.01074	25.50	No	
4	0.189	0.1187	0.207	0.2098	13.54	0.220	0.1458	0.00546	12.96	No	
5	0.220	0.1458	0.245	0.2654	17.13	0.267	0.1282	0.00991	23.53	No	
6	0.267	0.1282	0.280	0.1696	10.94	0.311	0.0000	0.00440	10.44	No	
7	0.983	0.0000	0.992	0.2560	16.52	0.999	0.0000	0.00217	5.15	No	

Track 15:

Type	Sample
Vial ID	MeOH blank
Description	MeOH Blank
Volume	2.0 µl

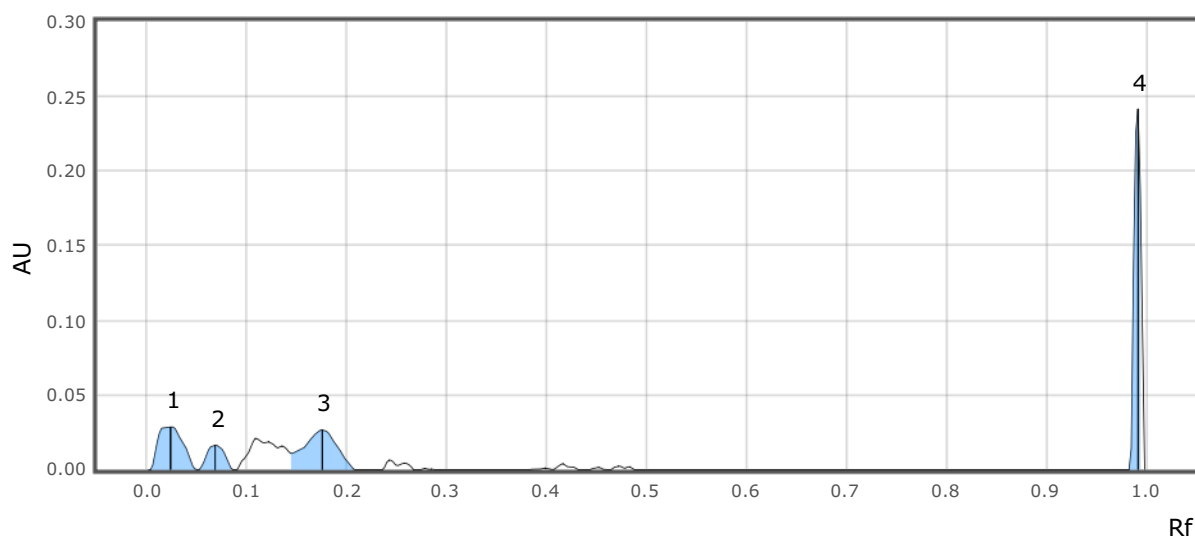


XHDa-re-3

visionCATS

Peak #	Start		Max			End		Area		Manual peak	Substance Name
	Rf	H	Rf	H	%	Rf	H	A	%		
1	0.004	0.0000	0.024	0.0285	9.12	0.050	0.0000	0.00082	19.70	No	
2	0.050	0.0000	0.068	0.0163	5.22	0.086	0.0000	0.00033	7.91	No	
3	0.144	0.0108	0.175	0.0266	8.51	0.209	0.0000	0.00102	24.39	No	
4	0.983	0.0000	0.992	0.2412	77.15	0.999	0.0000	0.00200	48.00	No	


Calibration results:

Height calibration for substance 8-THC @ RT White:

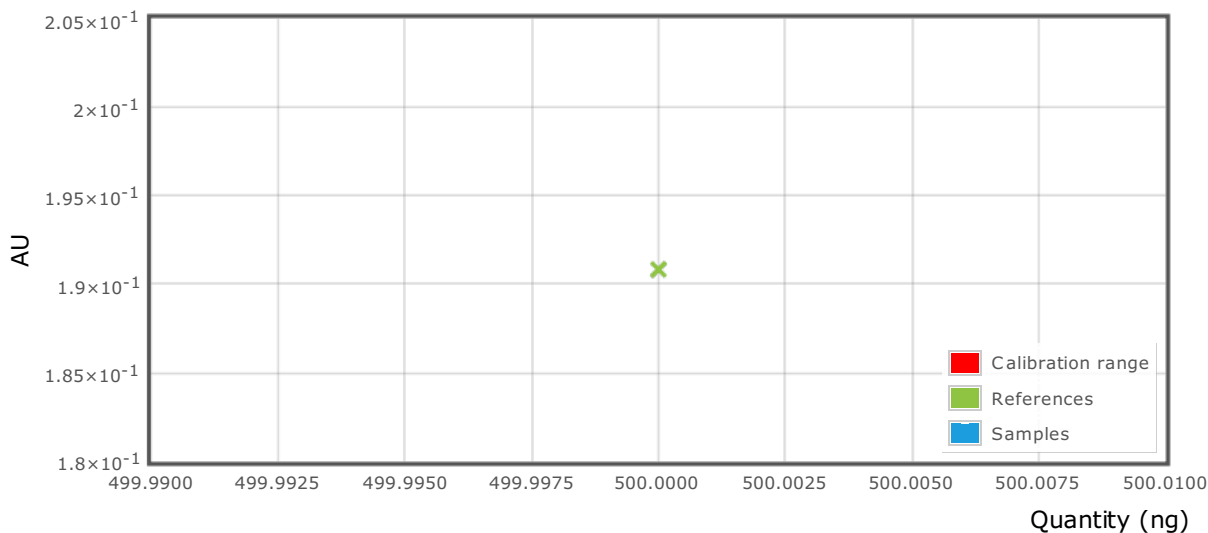


Regression mode	Linear-2
Range deviation	5.00 %
Related substances	Default
Number of references	1
Calibration function	$y=0x$
Coefficient of variation	CV 0.00 %
Correlation coefficient	n/a
	Unable to compute the results for this substance because there wasn't enough groups of references replicas (at least 1 for Linear-1, 2 for Linear2 and Mime-1 and 3 for Polynomial and MiMe-2)

Height calibration for substance 9-THC @ RT White:

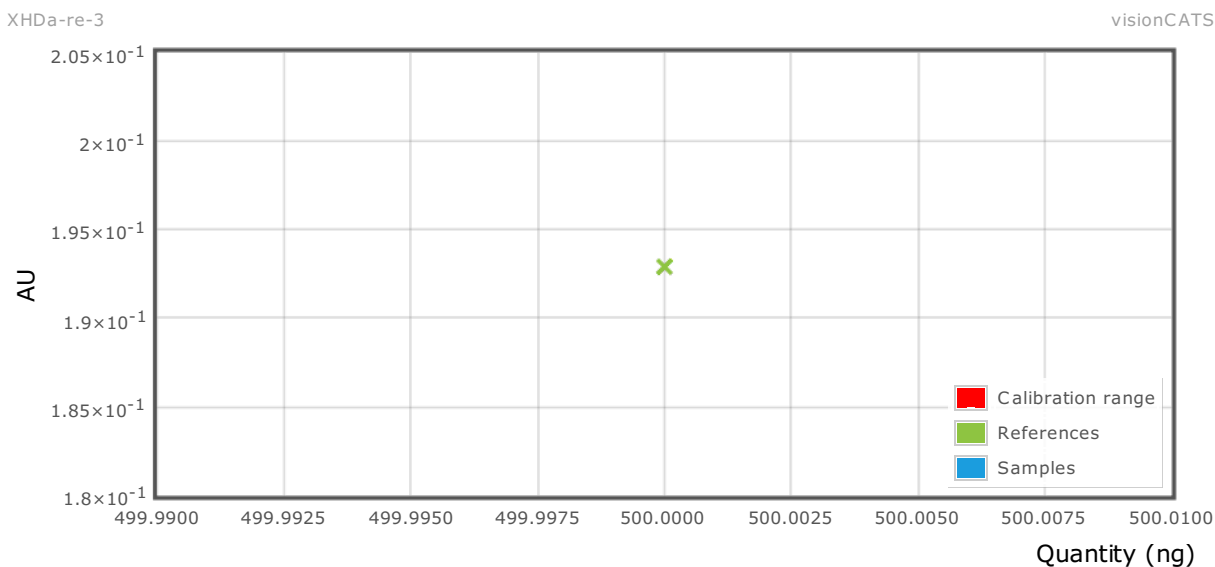

Regression mode	Linear-2
Range deviation	5.00 %
Related substances	Default
Number of references	1
Calibration function	$y=0x$
Coefficient of variation	CV 0.00 %
Correlation coefficient	n/a
	Unable to compute the results for this substance because there wasn't enough groups of references replicas (at least 1 for Linear-1, 2 for Linear2 and Mime-1 and 3 for Polynomial and MiMe-2)

Height calibration for substance CBC @ RT White:

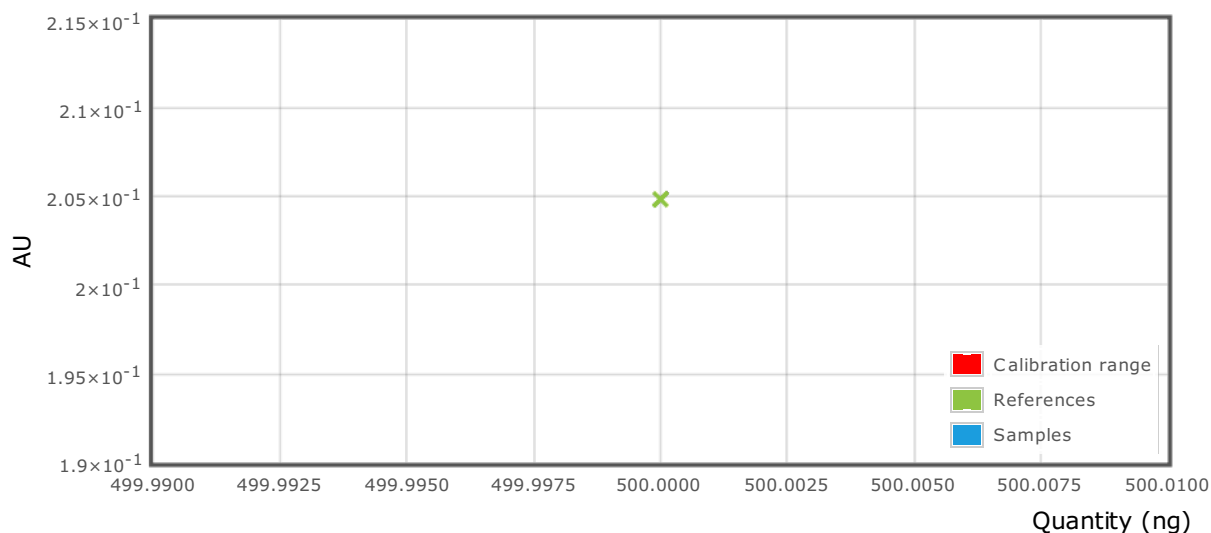


XHDa-re-3

visionCATS

Regression mode	Linear-2
Range deviation	5.00 %
Related substances	Default
Number of references	1
Calibration function	$y=0x$
Coefficient of variation	CV 0.00 %
Correlation coefficient	n/a
	Unable to compute the results for this substance because there wasn't enough groups of references replicas (at least 1 for Linear-1, 2 for Linear2 and Mime-1 and 3 for Polynomial and MiMe-2)

Height calibration for substance CBD @ RT White:

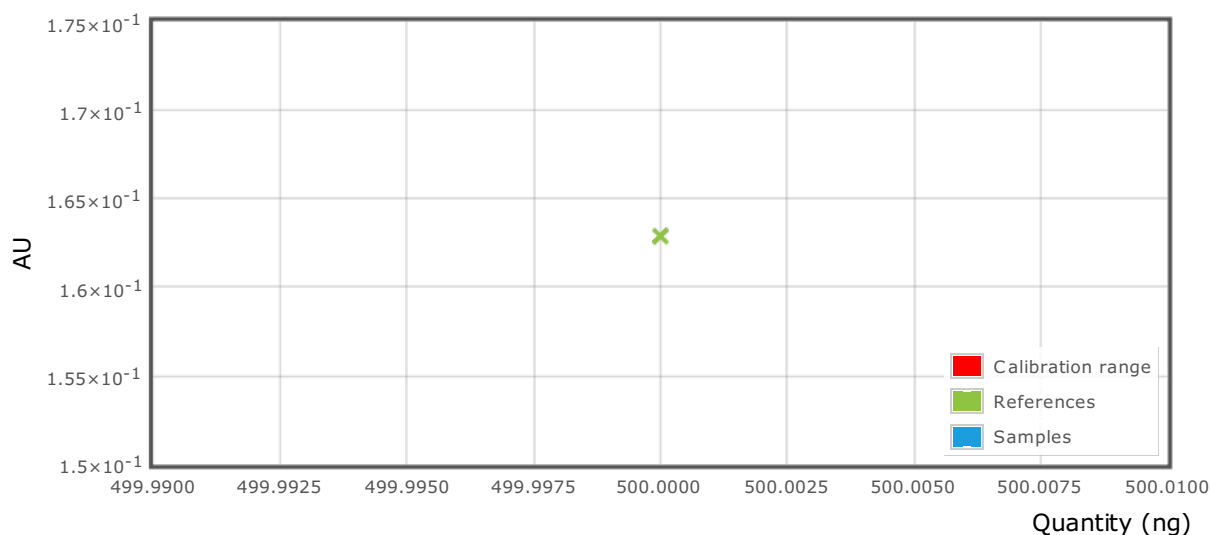


Regression mode	Linear-2
Range deviation	5.00 %
Related substances	Default
Number of references	1
Calibration function	$y=0x$
Coefficient of variation	CV 0.00 %
Correlation coefficient	n/a
	Unable to compute the results for this substance because there wasn't enough groups of references replicas (at least 1 for Linear-1, 2 for Linear2 and Mime-1 and 3 for Polynomial and MiMe-2)

Height calibration for substance CBDA @ RT White:


XHDa-re-3

visionCATS

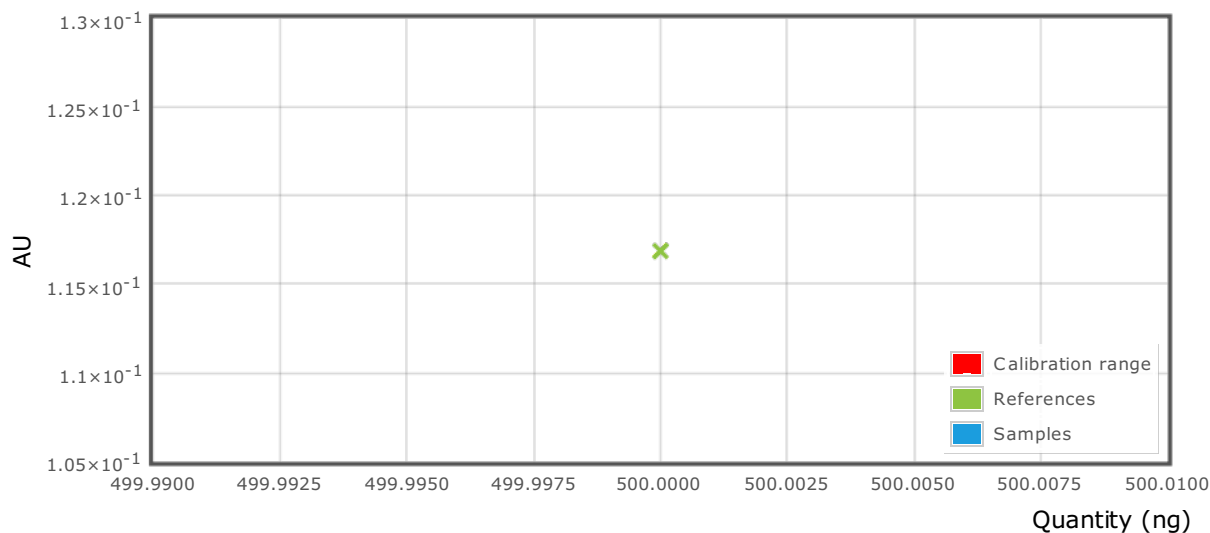

Regression mode	Linear-2
Range deviation	5.00 %
Related substances	Default
Number of references	1
Calibration function	$y=0x$
Coefficient of variation	CV 0.00 %
Correlation coefficient	n/a
	Unable to compute the results for this substance because there wasn't enough groups of references replicas (at least 1 for Linear-1, 2 for Linear2 and Mime-1 and 3 for Polynomial and MiMe-2)

Height calibration for substance CBDV @ RT White:

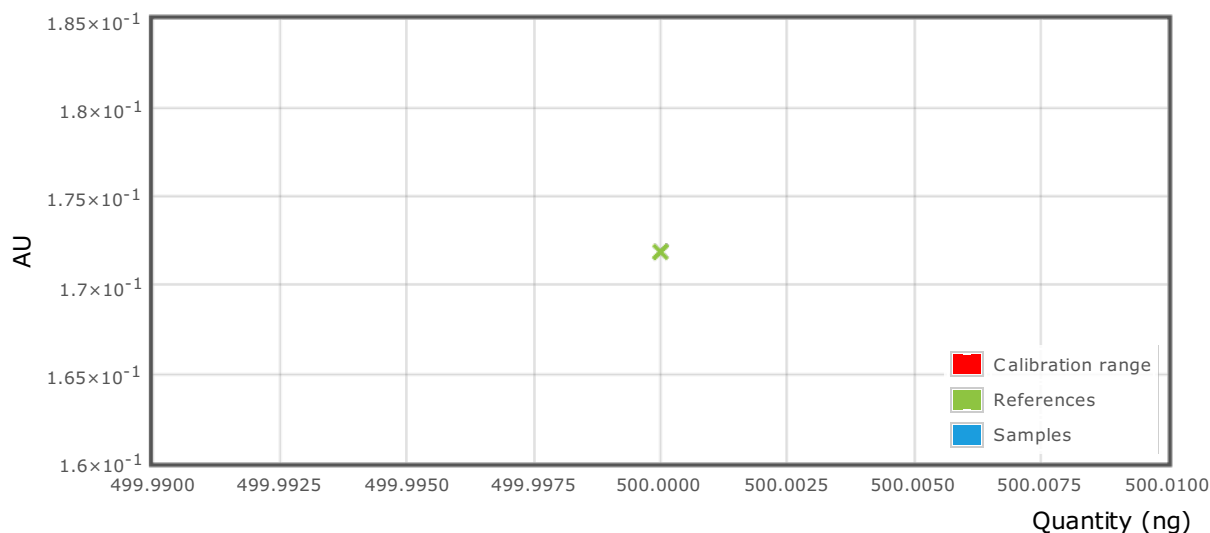


XHDa-re-3

visionCATS


Regression mode	Linear-2
Range deviation	5.00 %
Related substances	Default
Number of references	1
Calibration function	$y=0x$
Coefficient of variation	CV 0.00 %
Correlation coefficient	n/a
	Unable to compute the results for this substance because there wasn't enough groups of references replicas (at least 1 for Linear-1, 2 for Linear2 and Mime-1 and 3 for Polynomial and MiMe-2)

Height calibration for substance CBG @ RT White:

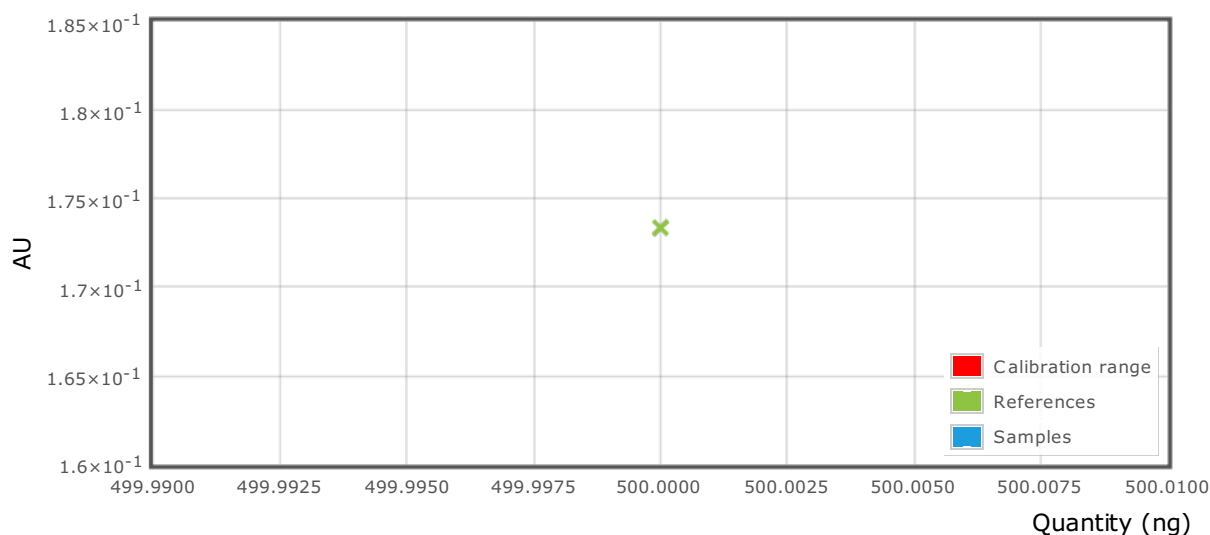


Regression mode	Linear-2
Range deviation	5.00 %
Related substances	Default
Number of references	1
Calibration function	$y=0x$
Coefficient of variation	CV 0.00 %
Correlation coefficient	n/a
	Unable to compute the results for this substance because there wasn't enough groups of references replicas (at least 1 for Linear-1, 2 for Linear2 and Mime-1 and 3 for Polynomial and MiMe-2)

Height calibration for substance CBGA @ RT White:

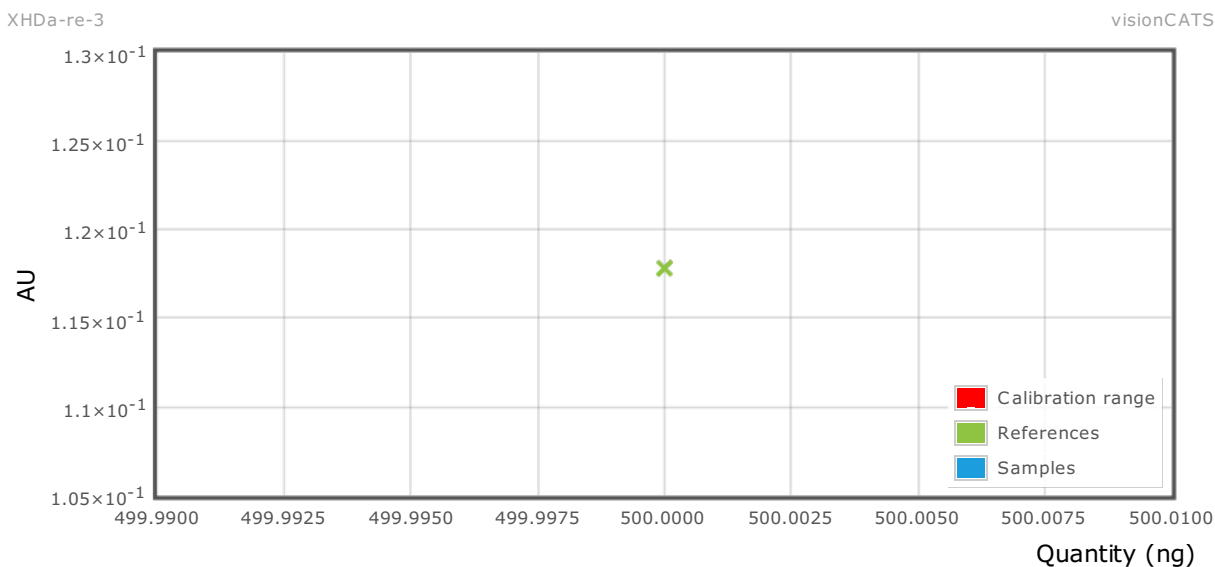

Regression mode	Linear-2
Range deviation	5.00 %
Related substances	Default
Number of references	1
Calibration function	$y=0x$
Coefficient of variation	CV 0.00 %
Correlation coefficient	n/a
	Unable to compute the results for this substance because there wasn't enough groups of references replicas (at least 1 for Linear-1, 2 for Linear2 and Mime-1 and 3 for Polynomial and MiMe-2)

Height calibration for substance CBN @ RT White:

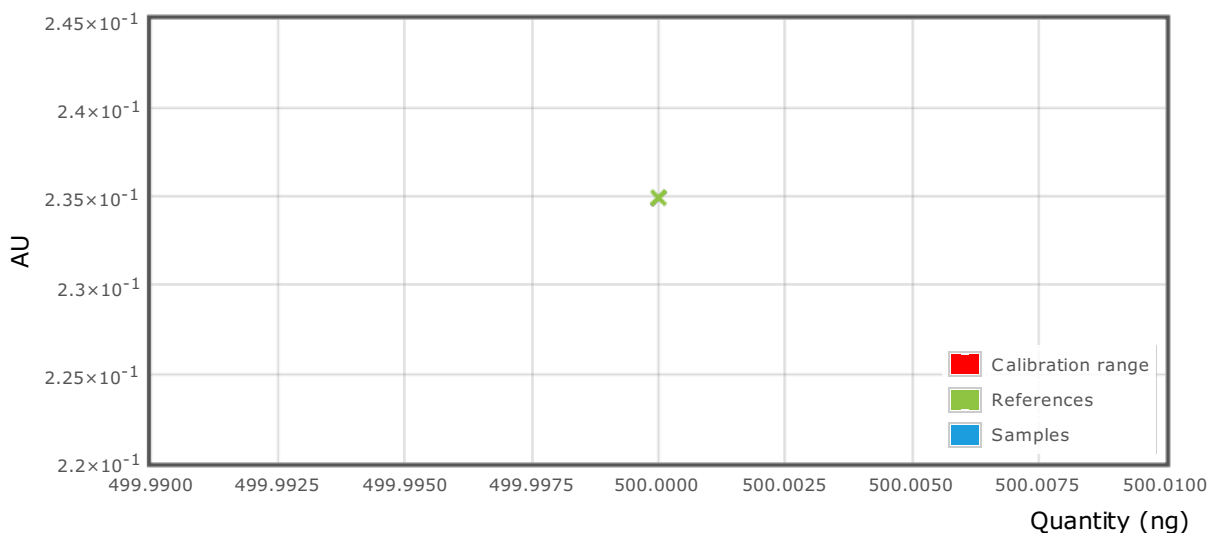


XHDa-re-3

visionCATS


Regression mode	Linear-2
Range deviation	5.00 %
Related substances	Default
Number of references	1
Calibration function	$y=0x$
Coefficient of variation	CV 0.00 %
Correlation coefficient	n/a
	Unable to compute the results for this substance because there wasn't enough groups of references replicas (at least 1 for Linear-1, 2 for Linear2 and Mime-1 and 3 for Polynomial and MiMe-2)

Height calibration for substance THCA-A @ RT White:

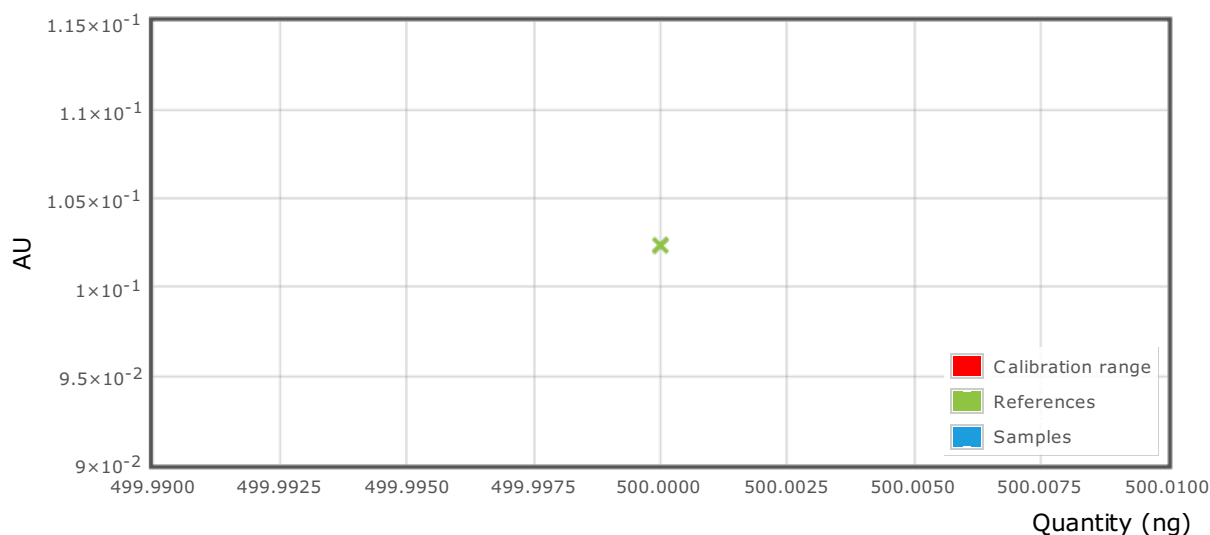

Regression mode	Linear-2
Range deviation	5.00 %
Related substances	Default
Number of references	1
Calibration function	$y=0x$
Coefficient of variation	CV 0.00 %
Correlation coefficient	n/a
	Unable to compute the results for this substance because there wasn't enough groups of references replicas (at least 1 for Linear-1, 2 for Linear2 and Mime-1 and 3 for Polynomial and MiMe-2)

Height calibration for substance THCV @ RT White:

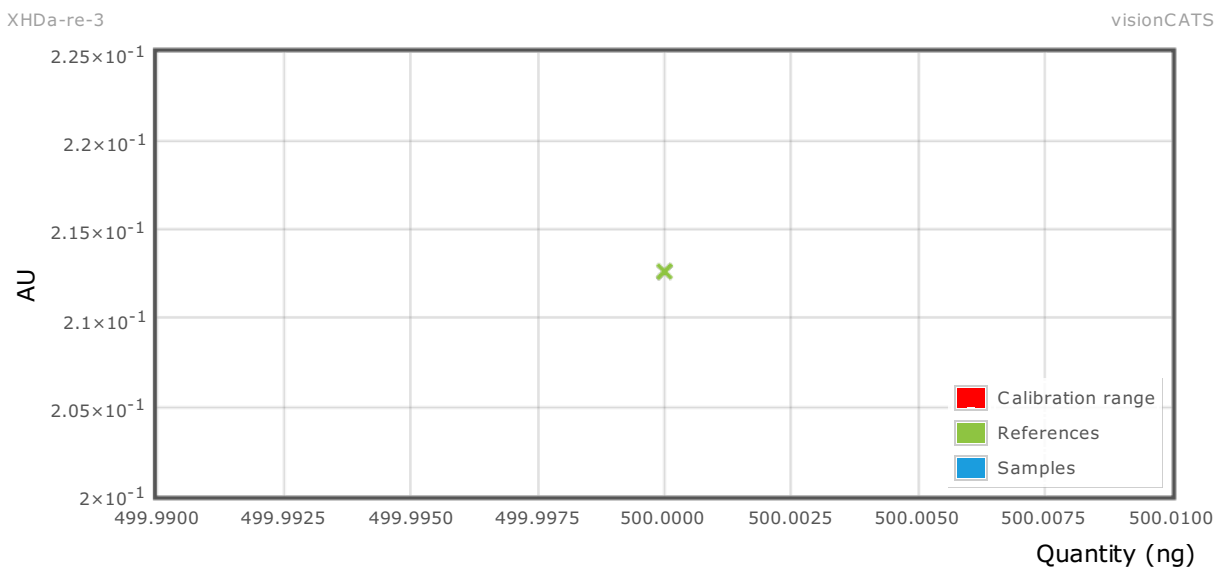

Regression mode	Linear-2
Range deviation	5.00 %
Related substances	Default
Number of references	1
Calibration function	$y=0x$
Coefficient of variation	CV 0.00 %
Correlation coefficient	n/a
	Unable to compute the results for this substance because there wasn't enough groups of references replicas (at least 1 for Linear-1, 2 for Linear2 and Mime-1 and 3 for Polynomial and MiMe-2)

Results:

Substance having no available results

	CBG	There wasn't any sample application available in the assignments for this substance. Please check that the peaks were correctly detected and assigned for this substance.
	CBDA	There wasn't any sample application available in the assignments for this substance. Please check that the peaks were correctly detected and assigned for this substance.
	CBD	There wasn't any sample application available in the assignments for this substance. Please check that the peaks were correctly detected and assigned for this substance.
	8-THC	There wasn't any sample application available in the assignments for this substance. Please check that the peaks were correctly detected and assigned for this substance.
	CBC	There wasn't any sample application available in the assignments for this substance. Please check that the peaks were correctly detected and assigned for this substance.
	CBGA	There wasn't any sample application available in the assignments for this substance. Please check that the peaks were correctly detected and assigned for this substance.
	THCV	There wasn't any sample application available in the assignments for this substance. Please check that the peaks were correctly detected and assigned for this substance.
	THCA-A	There wasn't any sample application available in the assignments for this substance. Please check that the peaks were correctly detected and assigned for this substance.
	CBN	There wasn't any sample application available in the assignments for this substance. Please check that the peaks were correctly detected and assigned for this substance.
	9-THC	There wasn't any sample application available in the assignments for this substance. Please check that the peaks were correctly detected and assigned for this substance.
	CBDV	There wasn't any sample application available in the assignments for this substance. Please check that the peaks were correctly detected and assigned for this substance.

A track marked with  means: this result is outside the regression range given by the reference assignments, but is included in the results because it is in the allowed range deviation.

Analyst:

Reviewer: